

GenieFluor 488 Anti-Mouse CD54 Antibody [YN1/1.7.4]
 Catalogue Code: AGEL0121

Antibody Data

Product SKU:	AGEL0121	Clone:	YN1/1.7.4
Applications:	FCM		
Reactivity:	Mouse		

Important Note:

Centrifuge before opening to ensure complete recovery of vial contents.

Product Information:

Alternate Names: Intercellular adhesion molecule 1;Icam1;MALA-2;MyD10;CD54;Icam-1;

Uniprot ID: P13597

Background: CD54 is a 90 kD immunoglobulin superfamily member also known as ICAM-1 and Ly-47. It is expressed on activated endothelial cells, high endothelial venules (HEV), T and B cells, monocytes/ macrophages, granulocytes, and dendritic cells. CD54 is an important intracellular adhesion molecule that participates in T cell-T cell, T cell-B cell, and T cell-target cell interactions via binding of LFA-1 (CD11a/CD18) and Mac-1 (CD11b/CD18). CD54 has also been shown to be involved in lymphocyte trafficking, making it an important molecule in many immune reactions and inflammation. CD54 is also a receptor for rhinovirus. The YN1/1.7.4 antibody has been reported to block binding of mouse CD54 to LFA-1 and Mac-1, inhibit cell-cell adhesion, and function in antigen presentation to T cells and leukocyte migration to inflammatory tissues.


Form: Liquid

Conjugation: Genie Fluor488

Size: 50 Tests, 100 Tests, 200 Tests

Host Species: Rat

Isotype: Rat IgG2b, κ

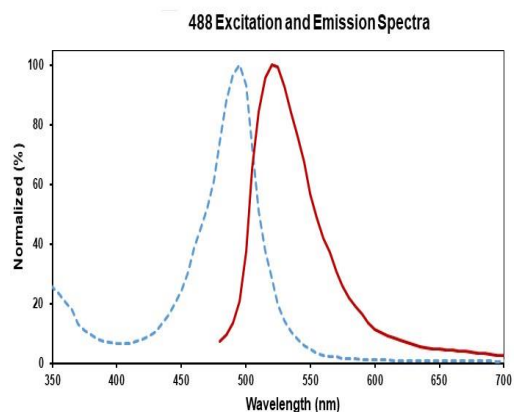

Isotype Control: Genie Fluor 488 Rat IgG2b, κ Isotype Control[LTF-2] [Product AGEL0121]

Storage Buffer: Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping: Biological ice pack at 4°C

Stability & Storage: Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

Recommended Usage: Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.