

**Antibody Data**

<b>Product SKU:</b>	<b>AGEL0674</b>	<b>Clone:</b>	<b>MIAP410</b>
<b>Applications:</b>	<b>FCM</b>		
<b>Reactivity:</b>	<b>Human;Mouse;Rat</b>		


**Important Note:**

Centrifuge before opening to ensure complete recovery of vial contents.

**Product Information:**

**Alternate Names:** Leukocyte surface antigen CD47;Cd47;Integrin-associated protein;IAP;  
**Uniprot ID:** Q08722 Q61735 P97829  
**Background:** CD47, also known as Integrin-Associated Protein (IAP), is a membrane protein of about 50 kD with an IgV-like extracellular domain, a five membrane-spanning segment and a short terminal cytoplasmic region. It is widely expressed on many cell types and often associated with beta 3 integrins. It has been reported that CD47 functions as a self marker. Red cells lacking CD47 were rapidly cleared from the bloodstream by splenic macrophages. By binding to SIRPα, CD47 controls hemostatic innate immune functions, such as phagocytosis and cell trafficking.

**Form:** Liquid  
**Conjugation:** PE/Cyanine 7  
**Size:** 20 Tests, 100 Tests, 200 Tests  
**Host Species:** Mouse  
**Isotype:** Mouse IgG1, κ

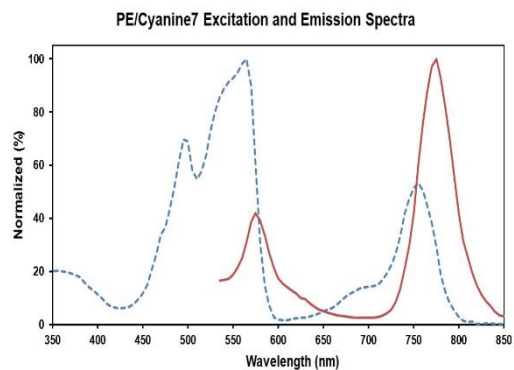

Ex:495;565;755 nm; Em:775 nm

**Isotype Control:** PE/Cyanine7 Mouse IgG1, κ Isotype Control[MOPC-21] [Product AGEL0674]  
**Storage Buffer:** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.  
**Shipping:** Biological ice pack at 4°C

**Stability & Storage:** Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

**Recommended Usage:** Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.