

Product Datasheet

PerCP/Cyanine5.5 Anti-Mouse CD127/IL-7RA Antibody [A7R34]

Catalogue Code: AGEL0743

Antibody Data

Product SKU: AGEL0743 Clone: A7R34

Applications: FCM

Reactivity: Mouse

Important Note:

Centrifuge before opening to ensure complete recovery of vial contents.

Product Information:

Alternate Names: Interleukin-7 receptor subunit alpha;IL-7 receptor subunit alpha;IL-7R subunit alpha;IL-7R subunit alpha;IL-7

7R-alpha;IL-7RA;CD127;

Uniprot ID: P16872

Background: CD127 is a 60-90 kD type I transmembrane glycoprotein also known as IL-7 receptor α

chain or IL-7Rα. It forms a heterodimer with the common γ chain (γc or CD132) which is shared with the receptors for IL-2, IL-4, IL-9, IL-13, IL-15, and IL-21. CD127 is expressed on immature B cells through early pre-B stage, thymocytes (except CD4/CD8 double positive thymocytes), peripheral T cells, and bone marrow stromal cells. CD127 has been reported to be an useful marker for identifying memory and effector T cells. The ligation of IL-7 with its receptor is important for stimulation of mature and immature T cells as well as

immature B cells proliferation and development.


Form: Liquid

Conjugation: PerCP/Cyanine 5.5

Size: 50 Tests, 100 Tests, 200 Tests

Host Species: Rat

Isotype: Rat IgG2a, κ

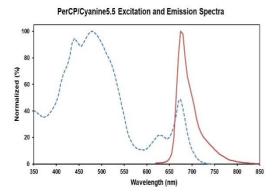

Ex:440;480;675 nm; Em:675 nm

Isotype Control: PerCP/Cyanine5.5 Rat IgG2a, κ Isotype Control[2A3] [Product AGEL0743]

Storage Buffer: Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping: Biological ice pack at 4°C



Stability & Storage:

Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

Recommended Usage:

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.