

**Antibody Data**

<b>Product SKU:</b>	<b>AGEL0775</b>	<b>Clone:</b>	<b>FGK4.5/FGK45</b>
<b>Applications:</b>	<b>FCM</b>		
<b>Reactivity:</b>	<b>Mouse</b>		

**Important Note:**


Centrifuge before opening to ensure complete recovery of vial contents.

**Product Information:**

**Alternate Names:** Tumor necrosis factor receptor superfamily member 5;Cd40;B-cell surface antigen CD40;Bp50;CD40L receptor;CD40;Tnfrsf5;  
**Uniprot ID:** P27512

**Background:** CD40 is a 48 kD type I transmembrane glycoprotein also known as Bp50. It is a member of the tumor necrosis factor receptor (TNFR) superfamily and is expressed on B cells, basal epithelial cells, macrophages, follicular dendritic cells, endothelial cells, and a subset of CD34+ hematopoietic progenitors. CD40 regulates B cell development/maturation, Ig isotype switching and, in combination with other signals such as IL-4, protects B cells from surface Ig-induced apoptosis and promotes proliferation. Interaction of CD40 with its ligand CD154 (gp39), which is expressed on activated T cells, is important in costimulation and immune regulation.

**Form:** Liquid  
**Conjugation:** FITC  
**Size:** 50 Tests, 100 Tests, 200 Tests  
**Host Species:** Rat  
**Isotype:** Rat IgG2a, κ

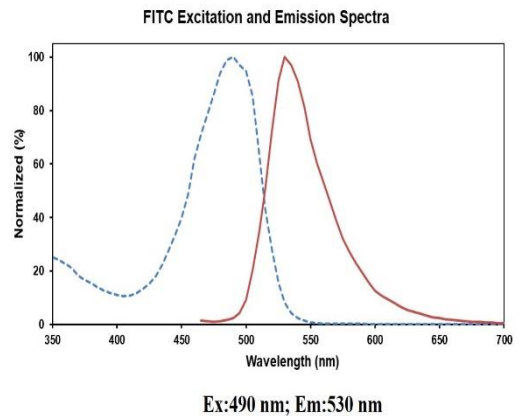

**Isotype Control:** FITC Rat IgG2a, κ Isotype Control[2A3] [Product AGEL0775]  
**Storage Buffer:** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.  
**Shipping:** Biological ice pack at 4°C

**Stability & Storage:** Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

**Recommended Usage:** Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.