

#### **Product Datasheet**

# PE/Cyanine7 Anti-Human CD45RA Antibody [HI100]

Catalogue Code: AGEL0958

## Antibody Data

Product SKU: AGEL0958 Clone: HI100

Applications: FCM

Reactivity: Human

# **Important Note:**

Centrifuge before opening to ensure complete recovery of vial contents.

### **Product Information:**

**Alternate Names:** Receptor-type tyrosine-protein phosphatase C;Ptprc;L-CA;Ly-5;T200;CD45;

Uniprot ID: P08575

**Background**: CD45RA is a 205-220 kD single chain type I glycoprotein. It is an exon 4 splice variant of

the tyrosine phosphatase CD45. The CD45RA isoform is expressed on resting/naïve T cells, medullary thymocytes, B cells and monocytes. CD45RA enhances both T cell receptor and B cell receptor signaling. CD45 non-covalently associates with lymphocyte phosphatase-associated phosphoprotein (LPAP) on T and B lymphocytes. CD45 has been reported to be associated with several other cell surface antigens including CD1, CD2, CD3, and CD4. CD45 has also been reported to bind galectin-1. CD45 isoform expression

can change in response to cytokines.


Form: Liquid

**Conjugation:** PE/Cyanine 7

Size: 20 Tests, 100 Tests, 200 Tests

Host Species: Mouse

**Isotype:** Mouse IgG2b, κ

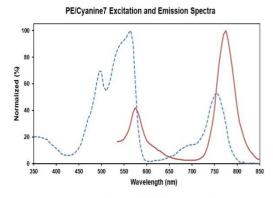

Ex:495;565;755 nm; Em:775 nm

**Isotype Control:** PE/Cyanine7 Mouse IgG2b, κ Isotype Control[MPC-11] [Product AGEL0958]

**Storage Buffer:** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

**Shipping:** Biological ice pack at 4°C



Stability & Storage:

Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

Recommended Usage:

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu L$  of antibody per test (million cells in 100  $\mu L$  staining volume or per 100  $\mu L$  of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.