

Product Datasheet

PE/Cyanine7 Anti-Human/Mouse CD44 Antibody [IM7]

Catalogue Code: AGEL1277

Antibody Data

Product SKU: AGEL1277 Clone: IM7

Applications: FCM

Reactivity: Human; Mouse

Important Note:

Centrifuge before opening to ensure complete recovery of vial contents.

Product Information:

Alternate Names: CD44 antigen;CD44;CDw44;Epican;Phagocytic glycoprotein 1;PGP-1;Phagocytic

glycoprotein I;PGP-I;CD44;LHR; MDU2; MDU3; MIC4;

Uniprot ID: P15379 P16070

Background: CD44 is a 80-95 kD glycoprotein also known as Hermes, Pgp1, H-CAM, or HUTCH. It is

expressed on all leukocytes, endothelial cells, hepatocytes, and mesenchymal cells. As B and T cells become activated or progress to the memory stage, CD44 expression increases from low or mid levels to high levels. Thus, CD44 has been reported to be a valuable marker for memory cell subsets. High CD44 expression on Treg cells has been associated with potent suppressive function via high production of IL-10. CD44 is an adhesion molecule involved in leukocyte attachment to and rolling on endothelial cells, homing to peripheral lymphoid organs and to the sites of inflammation, and leukocyte

aggregation.


Form: Liquid

Conjugation: PE/Cyanine 7

Size: 20 Tests, 100 Tests, 200 Tests

Host Species: Rat

Isotype: Rat IgG2b, κ

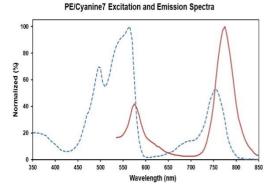

Ex:495;565;755 nm; Em:775 nm

Isotype Control: PE/Cyanine7 Rat IgG2b, κ Isotype Control[LTF-2] [Product AGEL1277]

Storage Buffer: Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping: Biological ice pack at 4°C



Stability & Storage: Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to

light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial

contents. This product is guaranteed up to one year from purchase.

Recommended Usage:

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.