

GenieFluor 647 Anti-Human/Mouse CD44 Antibody [IM7]

Catalogue Code: AGEL1280

Antibody Data

Product SKU:	AGEL1280	Clone:	IM7
Applications:	FCM		
Reactivity:	Human;Mouse		

Important Note:

Centrifuge before opening to ensure complete recovery of vial contents.


Product Information:

Alternate Names: CD44 antigen;CD44;CDw44;Epican;Phagocytic glycoprotein 1;PGP-1;Phagocytic glycoprotein I;PGP-I;CD44;LHR; MDU2; MDU3; MIC4; P15379 P16070

Uniprot ID:

Background: CD44 is a 80-95 kD glycoprotein also known as Hermes, Pgp1, H-CAM, or HUTCH. It is expressed on all leukocytes, endothelial cells, hepatocytes, and mesenchymal cells. As B and T cells become activated or progress to the memory stage, CD44 expression increases from low or mid levels to high levels. Thus, CD44 has been reported to be a valuable marker for memory cell subsets. High CD44 expression on Treg cells has been associated with potent suppressive function via high production of IL-10. CD44 is an adhesion molecule involved in leukocyte attachment to and rolling on endothelial cells, homing to peripheral lymphoid organs and to the sites of inflammation, and leukocyte aggregation.

Form:	Liquid
Conjugation:	Genie Fluor647
Size:	20 Tests, 100 Tests, 200 Tests
Host Species:	Rat
Isotype:	Rat IgG2b, κ

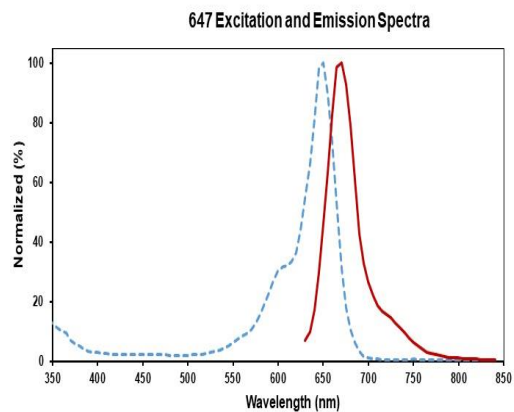

Isotype Control:	Genie Fluor 647 Rat IgG2b, κ Isotype Control[LTF-2] [Product AGEL1280]
Storage Buffer:	Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.
Shipping:	Biological ice pack at 4°C

Stability & Storage: Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

Recommended Usage: Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.
