

Product Datasheet

PE Anti-Human CD45RO Antibody [UCHL1]

Catalogue Code: AGEL1697

Antibody Data

Product SKU: AGEL1697 Clone: UCHL1

Applications: FCM

Reactivity: Human

Important Note:

Centrifuge before opening to ensure complete recovery of vial contents.

Product Information:

Alternate Names: Uniprot ID: -

Background: CD45RO is a 180 kD single chain membrane glycoprotein. It is a splice variant of tyrosine

phosphatase CD45, lacking the A, B, and C determinants. The CD45RO isoform is expressed on activated and memory T cells, some B cell subsets, activated monocytes/macrophages, and granulocytes. CD45RO enhances both T cell receptor and B cell receptor signaling mediated activation. CD45 and its isoforms non-covalently associate with lymphocyte phosphatase-associated phosphoprotein (LPAP) on T and B lymphocytes. CD45 has been reported to be associated with several other cell surface antigens including CD1, CD2, CD3, and CD4. CD45 has also been reported to bind galectin-1 and CD22. CD45 isoform expression can change in response to cytokines.


Form: Liquid

Conjugation: PE

Size: 20 Tests, 100 Tests, 200 Tests

Host Species: Mouse

Isotype: Mouse IgG2a, κ

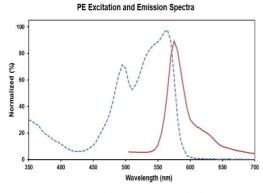

Ex:495;565 nm; Em:575 nm

Isotype Control: PE Mouse IgG2a, κ Isotype Control[C1.18.4] [Product AGEL1697]

Storage Buffer: Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping: Biological ice pack at 4°C



Stability & Storage:

Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

Recommended Usage:

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.