

Antibody Data

Product SKU:	AGEL2037	Clone:	TC11-18H10.1
Applications:	ICFCM		
Reactivity:	Mouse		

Important Note:

Centrifuge before opening to ensure complete recovery of vial contents.

Product Information:

Alternate Names: Interleukin-17A;IL-17;IL-17A;Cytotoxic T-lymphocyte-associated antigen 8;CTLA-8;CTLA8;

Uniprot ID: Q62386

Background: IL-17, also known as CTLA-8, is a T cell-expressed pleiotropic cytokine that exhibits a high degree of homology to a protein encoded by the ORF13 gene of herpes virus Saimiri. IL-17 is produced by Th cells (Th17) that are distinct from the traditional Th1- and Th2-cell subsets. IL-23 plays an important role in triggering IL-17 production. Both recombinant and natural IL-17 have been shown to exist as disulfide linked homodimers. IL-17 exhibits multiple biological activities on a variety of cells including: the induction of IL-6 and IL-8 production in fibroblasts, activation of NF-κB, and costimulation of T cell proliferation. IL-17 is an essential inflammatory mediator in the development of autoimmune diseases. Neutralization of IL-17 with monoclonal antibody is able to ameliorate the disease course.


Form: Liquid

Conjugation: APC

Size: 50 Tests, 100 Tests, 200 Tests

Host Species: Rat

Isotype: Rat IgG1, κ

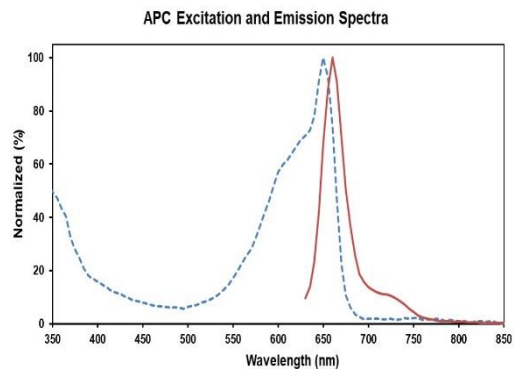

Ex:650 nm; Em:660 nm

Isotype Control: APC Rat IgG1, κ Isotype Control[HRPN] [Product AGEL2037]

Storage Buffer: Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping: Biological ice pack at 4°C

Stability & Storage: Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

Recommended Usage: Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

