

Antibody Data

Product SKU:	AGEL2277	Clone:	A019D5
Applications:	FCM		
Reactivity:	Human		

Important Note:

Centrifuge before opening to ensure complete recovery of vial contents.

Product Information:

Alternate Names: Interleukin-7 receptor subunit alpha;IL-7RA;CDw127;CD127;IL-7R α ;

Uniprot ID: P16871

Background: CD127 is a 60-90 kD type I transmembrane glycoprotein also known as IL-7 receptor α chain or IL-7R α . It forms a heterodimer with the common γ chain (γ c or CD132) which is shared with the receptors for IL-2, IL-4, IL-9, IL-13, IL-15, and IL-21. CD127 is expressed on immature B cells through early pre-B stage cells, thymocytes (except CD4/CD8 double positive thymocytes), peripheral T cells, and bone marrow stromal cells. CD127 has been reported to be a useful marker for identifying memory and effector T cells. Studies have shown that CD127 expression is down-modulated on Treg cells. It can be used as a marker for differentiation of Treg and conventional T cells. The ligation of IL-7 with its receptor is important for stimulation of mature and immature T cells as well as immature B cell proliferation and development.


Form: Liquid

Conjugation: PerCP/Cyanine 5.5

Size: 20 Tests, 100 Tests, 200 Tests

Host Species: Mouse

Isotype: Mouse IgG1, κ

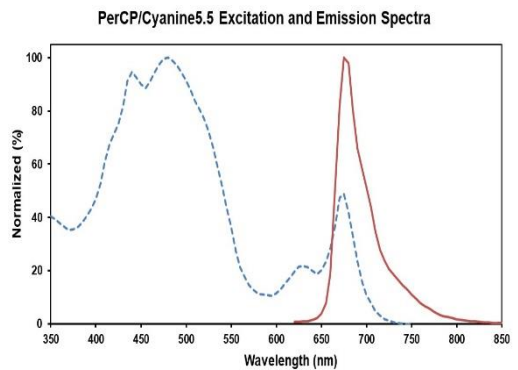

Ex:440;480;675 nm; Em:675 nm

Isotype Control: PerCP/Cyanine5.5 Mouse IgG1, κ Isotype Control[MOPC-21] [Product AGEL2277]

Storage Buffer: Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping: Biological ice pack at 4°C

Stability & Storage: Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

Recommended Usage: Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.