

**GenieFluor 647 Anti-Human CD20 Antibody [2H7]**  
 Catalogue Code: AGEL2400

**Antibody Data**

<b>Product SKU:</b>	<b>AGEL2400</b>	<b>Clone:</b>	<b>2H7</b>
<b>Applications:</b>	<b>FCM</b>		
<b>Reactivity:</b>	<b>Human</b>		

**Important Note:**

Centrifuge before opening to ensure complete recovery of vial contents.

**Product Information:**

**Alternate Names:** Bp35;Leukocyte surface antigen Leu-16;MS4A1;B-lymphocyte surface antigen B1;

**Uniprot ID:** P11836

**Background:** CD20 is a 33-37 kD, four transmembrane spanning protein, also known as B1 and Bp35. CD20 is expressed on pre-B-cells, resting and activated B cells (not plasma cells), some follicular dendritic cells, and at low levels on a T cell subset. CD20 is heavily phosphorylated on activated B cells and malignant B cells. Homo-oligomeric complexes of CD20 are thought to form Ca<sup>2+</sup> conductive ion channels in the plasma membrane of B cells. The CD20 molecule is involved in B-cell activation and is associated with various Src family kinases (Lyn, Lck, Fyn). It exists in a complex with MHC class I and II, CD53, CD81, and CD82.


**Form:** Liquid

**Conjugation:** Genie Fluor647

**Size:** 20 Tests, 100 Tests, 200 Tests

**Host Species:** Mouse

**Isotype:** Mouse IgG2b, κ

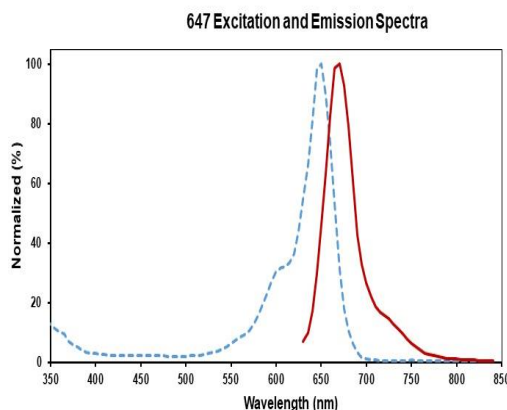

**Isotype Control:** Genie Fluor 647 Mouse IgG2b, κ Isotype Control[MPC-11] [Product AGEL2400]

**Storage Buffer:** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

**Shipping:** Biological ice pack at 4°C

**Stability & Storage:** Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

**Recommended Usage:** Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.