

#### **Product Datasheet**

# PE Anti-Human CD133 Antibody [W6B3C1]

Catalogue Code: AGEL2730

## Antibody Data

Product SKU: AGEL2730 Clone: W6B3C1

Applications: FCM

Reactivity: Human

#### **Important Note:**

Centrifuge before opening to ensure complete recovery of vial contents.

## **Product Information:**

**Alternate Names:** Prominin-1;PROM1;AC133;PROML1;

Uniprot ID: 043490

**Background**: CD133, also known as Prominin-1 and AC133 antigen, is a 120 kD pentaspan glycoprotein

with 5 transmembrane domains. CD133 was initially described as a surface antigen specific for human hematopoietic stem cells and as a marker for murine neuroepithelial cells and some embryonic epithelia. Later on, CD133 was found on other stem cells, including endothelial progenitor cells, glioblastomas, neuronal, and glial stem cells. In addition to stem cells for normal tissue, CD133 was found on cancer cells, such as some leukemia cells and brain tumor cells. Although the biological function of CD133 is not completely understood, CD133 has been extensively used as a stem cell marker for

normal and cancerous tissues.


Form: Liquid

Conjugation: PE

Size: 20 Tests, 100 Tests, 200 Tests

Host Species: Mouse

**Isotype:** Mouse IgG1, κ

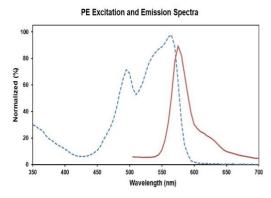

Ex:495;565 nm; Em:575 nm

**Isotype Control:** PE Mouse IgG1, κ Isotype Control[MOPC-21] [Product AGEL2730]

**Storage Buffer:** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

**Shipping:** Biological ice pack at 4°C



**Stability & Storage:** Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to

light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial

contents. This product is guaranteed up to one year from purchase.

Recommended Usage:

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.