

**GenieFluor Red 780 Anti-Human CD11b Antibody [ICRF44]**  
 Catalogue Code: AGEL2997

**Antibody Data**

<b>Product SKU:</b>	<b>AGEL2997</b>	<b>Clone:</b>	<b>ICRF44</b>
<b>Applications:</b>	<b>FCM</b>		
<b>Reactivity:</b>	<b>Human</b>		


**Important Note:**

Centrifuge before opening to ensure complete recovery of vial contents.

**Product Information:**

**Alternate Names:** Integrin alpha-M;Itgam;CD11 antigen-like family member B;CR-3 alpha chain;Leukocyte adhesion receptor MO1;CD11b;  
**Uniprot ID:** P11215  
**Background:** CD11b is a 165-170 kD type I transmembrane glycoprotein also known as  $\alpha$ MIntegrin, Mac-1, CR3, and C3biR. CD11b non-covalently associates with integrin  $\beta$ 2(CD18) and is expressed on granulocytes, monocytes/macrophages, dendritic cells, NK cells, and subsets of T and B cells. CD11b/CD18 is critical for the transendothelial migration of monocytes and neutrophils. It is also involved in granulocyte adhesion, phagocytosis, and neutrophil activation. CD11b/CD18 interacts with ICAM-1 (CD54), ICAM-2 (CD102), ICAM-4, CD14, CD23, heparin, iC3b, fibrinogen, and factor X.

**Form:** Liquid  
**Conjugation:** Genie FluorRed 780  
**Size:** 20 Tests, 100 Tests, 200 Tests  
**Host Species:** Mouse  
**Isotype:** Mouse IgG1,  $\kappa$

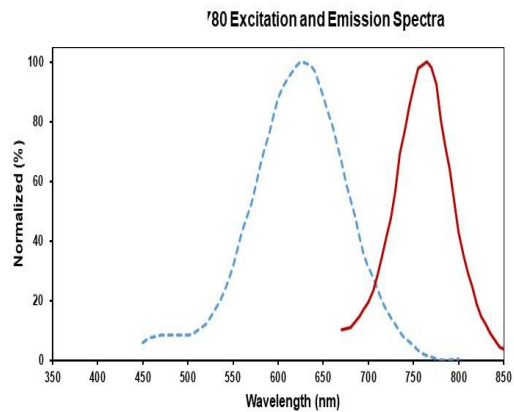

**Isotype Control:** Genie Fluor Red 780 Mouse IgG1,  $\kappa$  Isotype Control[MOPC-21] [Product AGEL2997]  
**Storage Buffer:** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.  
**Shipping:** Biological ice pack at 4°C

**Stability & Storage:** Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.

**Recommended Usage:** Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.