

Product Datasheet **FITC Anti-Human CD56 Antibody** [5.1H11] Catalogue Code: AGEL3041

## Antibody Data


Product SKU:	AGEL3041	Clone:	5.1H11	
Applications:	FCM			
Reactivity:	Human			

## Important Note:

Centrifuge before opening to ensure complete recovery of vial contents.

## **Product Information:**

Alternate Names: Uniprot ID:	Neural cell adhesion molecule 1;NCAM1;N-CAM-1;NCAM-1;NCAM; P13591		
Background:	CD56 is a single transmembrane glycoprotein also known as NCAM (Neural Cell Adhesion Molecule), Leu-19, or NKH1. It is a member of the Ig superfamily. The 140 kD isoform is expressed on NK cells and NK-T cells. CD56 is also expressed in the brain (cerebellum and cortex) and at neuromuscular junctions. Certain large granular lymphocyte (LGL) leukemias, small-cell lung carcinomas, neuronal derived tumors, myelomas, and myeloid leukemias also express CD56. CD56 plays a role in homophilic and heterophilic adhesion via binding to itself or heparin sulfate.		
Form:	Liquid	FITC Excitation and Emission Spectra	
Conjugation:	FITC	100 -	
Size:	20 Tests, 100 Tests, 200 Tests		
Host Species:	Mouse	(%) 40	
Isotype: Mo	Mouse IgG1, κ	2 2 0 3 5 4 0 4 5 5 5 5 5 6 6 0 6 5 7 0 Wavelength (nm)	

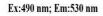

**Isotype Control:** FITC Mouse IgG1, κ Isotype Control[MOPC-21] [Product AGEL3041]

**Storage Buffer:** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping: Biological ice pack at 4°C



- **Stability & Storage:** Keep as concentrated solution. Store at 2~8°C and protected from prolonged exposure to light. Do not freeze. Centrifuge before opening to ensure complete recovery of vial contents. This product is guaranteed up to one year from purchase.
- **Recommended** Usage: Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.