MMP3 Rabbit Polyclonal Antibody



CAB1202

Product Information

Size:

20uL, 50uL, 100uL, 200uL

Observed MW:

54kDa

Calculated MW:

53kDa

Applications:

WB IF

Reactivity:

Human, Mouse, Rat

Protein Background

Proteins of the matrix metalloproteinase (MMP) family are involved in the breakdown of extracellular matrix in normal physiological processes, such as embryonic development, reproduction, and tissue remodeling, as well as in disease processes, such as arthritis and metastasis. Most MMP's are secreted as inactive proproteins which are activated when cleaved by extracellular proteinases. This gene encodes an enzyme which degrades fibronectin, laminin, collagens III, IV, IX, and X, and cartilage proteoglycans. The enzyme is thought to be involved in wound repair, progression of atherosclerosis, and tumor initiation. The gene is part of a cluster of MMP genes which localize to chromosome 11q22.3.

Immunogen information

Gene ID:

4314

Uniprot

P08254

Synonyms:

MMP3; CHDS6; MMP-3; SL-1; STMY; STMY1; STR1 **Antibody Information**

Recommended dilutions:

WB 1:500 - 1:2000 IF 1:50 -

1:100

Source:

Rabbit

Isotype:

Immunogen:

Recombinant fusion protein containing a sequence corresponding to amino acids 248-477 of human MMP3 (NP_002413.1).

Storage:


Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.02% IgG

sodium azide, 50% glycerol, pH7.3.


Purification:

Affinity purification


Product Images

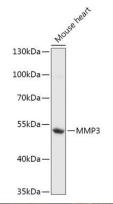

Western blot analysis of extracts of Mouse heart, using MMP3 antibody (CAB1202) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (CABS014) at 1:10000 dilution. Lysates/proteins: 25ug per lane. Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (CABM00020). Exposure time: 30s.

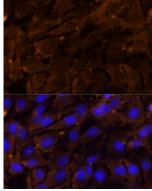

Immunofluorescence analysis of C6 cells using MMP3 Polyclonal Antibody (CAB1202) at dilution of 1:100 (40x lens). Blue: DAPI for nuclear staining.

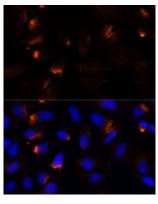

Immunofluorescence analysis of U-2 OS cells using MMP3 Polyclonal Antibody (CAB1202) at dilution of 1:100 (40x lens). Blue: DAPI for nuclear staining.