

# Human Trail R4/DcR2 PharmaGenie ELISA Kit

SKU: HUDC0020

**Technical Manual** 

For research use only



# **Table of Contents**

1.	Intended use	3
2.	Introduction	3
2.1.	Summary	3
2.2.	Principle of the method	3
3.	Reagents provided and reconstitution	4
4.	Materials required but not provided	4
5.	Storage Instructions	4
6.	Specimen collection, processing & storage	5
7.	Safety & precautions for use	5
8.	Assay Preparation	6
8.1.	Assay Design	6
8.2.	Preparation of Wash Buffer	6
8.3.	Preparation of Standard Diluent Buffer	6
8.4.	Preparation of Standard	7
8.5.	Preparation of Samples	7
8.6.	Preparation of Biotinylated Anti TRAIL R4	7
8.7.	Preparation of HRP-Conjugate	8
9.	Method	9
10.	Data Analysis	10
11.	Assay limitations	10
12.	Performance Characteristics	11
12.	1. Sensitivity	11
12.2	2. Specificity	11
12.3	3. Precision	11
12.4	4. Dilution Parallelism	11
12.5	5. Spike Recovery	11
12.6	6. Expected Serum Values	11
13.	Assay Summary	12



# **Human CD264 PharmaGenie ELISA Kit**

#### 1. Intended use

The ELISA Genie Human CD264 PharmaGenie ELISA Kit is a solid phase sandwich ELISA for the *in-vitro* qualitative and quantitative determination of TRAIL R4/DcR2 in cell culture supernatants, buffered solutions or human serum, plasma, or other body fluids. This assay will recognize both natural and recombinant human TRAIL R4/DcR2.

PharmaGenie ELISA Kits from ELISA Genie are a premium range of pre-coated ELISA kits especially designed for scientists working in pharmaceutical, biotech & CRO sectors. PharmaGenie ELISA kits are produced using high quality monoclonal antibody pairs & optimized reagents that have been manufactured according ISO 9001:2000 quality systems and are excellent assays to help discover our future.

This kit has been configured for research use only.

#### 2. Introduction

## 2.1. Summary

Human Trail R4, also called DcR1, CD264 or TRUNND is a type 1 TNF R family membrane protein which is a receptor for TRAIL. In the Trail receptor family, Trail R4 (DcR2) but also Trail R3 (DcR1) antagonize TRAIL-induced apoptosis whereas Trail R1 (DR4) and Trail R2 (DR5) transducer an apoptosis signal. Trail R4 has a truncated intracellular domain. Overexpression of Trail R4 has been shown to protect cells bearing Trail R1 and/or Trail R2 from Trail-induced apoptosis. Many tumor cell lines don't express the decoy receptors and are therefore sensitive to TRAIL if they express Trail R1 and/or Trail R2.

#### 2.2. Principle of the method

A capture Antibody highly specific for TRAIL R4 has been coated to the wells of the microtitre strip plate provided during manufacture. Binding of TRAIL R4 in samples and known standards to the capture antibodies is completed and then any excess unbound analyte is removed. During the next incubation period the binding of the biotinylated anti- TRAIL R4 secondary antibody to the analyte occurs. Any excess unbound secondary antibody is then removed. The HRP conjugate solution is then added to every well including the zero wells, following incubation excess conjugate is removed by careful washing. A chromogen substrate is added to the wells resulting in the progressive development of a blue coloured complex with the conjugate. The colour development is then stopped by the addition of acid turning the resultant final product yellow. The intensity of the produced coloured complex is directly proportional to the concentration of TRAIL R4 present in the samples and standards. The absorbance of the colour complex is then measured and the generated OD values for each standard are plotted against expected concentration forming a standard curve. This standard curve can then be used to accurately determine the concentration of TRAIL R4 in any sample tested.



# 3. Reagents provided and reconstitution

Reagents (Store@2- 8°C)	Quantity 1x48 well kit	Quantity 1x96 well kit	Quantity 2x96 well kit	Reconstitution
96 well microtiter strip plate	1/2	1	2	Ready to use (Pre-coated)
Plastic plate covers	2	2	4	n/a
Standard: 10 000pg/ml	1 vial	2 vials	4 vials	Reconstitute as directed on the vial (see Assay preparation, section 8)
Standard Diluent (Buffer)	1 vial (25ml)	1 vial (25ml)	1 vial (25ml)	10x Concentrate, dilute in distilled water (see reagent preparation, section 8)
Biotinylated anti- TRAIL R4	1 vial (0.4ml)	1 vial (0.4ml)	2 vials (0.4ml)	Dilute in Biotinylated Antibody Diluent (see Assay preparation, section 8)
Biotinylated Antibody Diluent	1 vial (7ml)	1 vial (7ml)	1 vial (13ml)	Ready to use
Streptavidin-HRP	1 vial (5µl)	2 vials (5µl)	4 vials (5µl)	Add 0.5ml of HRP diluent prior to use (see Assay preparation, section 8)
HRP Diluent	1 vial (23ml)	1 vial (23ml)	1 vial (23ml)	Ready to use
Wash Buffer	1 vial (10ml)	1 vial (10ml)	2 vials (10ml)	200x Concentrate dilute in distilled water (see Assay preparation, section 8)
TMB Substrate	1 vial (11ml)	1 vial (11ml)	1 vial (24ml)	Ready to use
H <sub>2</sub> SO <sub>4</sub> stop reagent	1 vial (11ml)	1 vial (11ml)	2 vials (11ml)	Ready to use

# 4. Materials required but not provided

- Microtiter plate reader fitted with appropriate filters (450nm required with optional 630nm reference filter)
- Microplate washer or wash bottle
- 10, 50, 100, 200 and 1,000µl adjustable single channel micropipettes with disposable tips
- 50-300µl multi-channel micropipette with disposable tips
- Multichannel micropipette reagent reservoirs
- · Distilled water
- Vortex mixer
- Miscellaneous laboratory plastic and/or glass, if possible sterile

# 5. Storage Instructions

Store kit reagents between 2 and 8°C. Immediately after use remaining reagents should be returned to cold storage (2-8°C). Expiry of the kit and reagents is stated on box front labels. The expiry of the kit components can only be guaranteed if the components are stored properly, and if, in case of repeated use of one component, the reagent is not contaminated by the first handling.

Wash Buffer: Once prepared store at 2-8°C for up to 1 week

Standard Diluent Buffer: Once prepared store at 2-8°C for up to 1 week

Standards: Once prepared use immediately and do not store

Biotinylated Secondary Antibody: Once prepared use immediately and do not store

Streptavidin-HRP: Once prepared use immediately and do not store



# 6. Specimen collection, processing & storage

Cell culture supernatants, serum, plasma or other biological samples will be suitable for use in the assay. Remove serum from the clot or red cells, respectively, as soon as possible after clotting and separation.

**Cell culture supernatants**: Remove particulates and aggregates by spinning at approximately 1000 x g for 10 min.

**Serum**: Use pyrogen/endotoxin free collecting tubes. Serum should be removed rapidly and carefully from the red cells after clotting. Following clotting, centrifuge at approximately 1000 x g for 10 min and remove serum.

Plasma: EDTA, citrate and heparin plasma can be assayed. Spin samples at 1000 x g for 30 min to remove

**Storage**: If not analyzed shortly after collection, samples should be aliquoted (250-500µI) to avoid repeated freeze-thaw cycles and stored frozen at -70°C. Avoid multiple freeze-thaw cycles of frozen specimens.

**Recommendation:** Do not thaw by heating at 37°C or 56°C. Thaw at room temperature and make sure that sample is completely thawed and homogeneous before use. When possible avoid use of badly haemolysed or lipemic sera. If large amounts of particles are present these should be removed prior to use by centrifugation or filtration.

## 7. Safety & precautions for use

- Handling of reagents, serum or plasma specimens should be in accordance with local safety procedures,
   e.g.CDC/NIH Health manual: "Biosafety in Microbiological and Biomedical Laboratories" 1984
- Laboratory gloves should be worn at all times
- Avoid any skin contact with H<sub>2</sub>SO<sub>4</sub> and TMB. In case of contact, wash thoroughly with water
- · Do not eat, drink, smoke or apply cosmetics where kit reagents are used
- Do not pipette by mouth
- When not in use, kit components should be stored refrigerated or frozen as indicated on vials or bottles labels
- All reagents should be warmed to room temperature before use. Lyophilized standards should be discarded after use
- Once the desired number of strips has been removed, immediately reseal the bag to protect the remaining strips from deterioration
- · Cover or cap all reagents when not in use
- Do not mix or interchange reagents between different lots
- Do not use reagents beyond the expiration date of the kit
- Use a clean disposable plastic pipette tip for each reagent, standard, or specimen addition in order to avoid cross contamination, for the dispensing of H<sub>2</sub>SO<sub>4</sub> and substrate solution, avoid pipettes with metal parts
- Use a clean plastic container to prepare the washing solution
- Thoroughly mix the reagents and samples before use by agitation or swirling
- All residual washing liquid must be drained from the wells by efficient aspiration or by decantation followed by tapping the plate forcefully on absorbent paper. Never insert absorbent paper directly into the wells
- The TMB solution is light sensitive. Avoid prolonged exposure to light. Also, avoid contact of the TMB solution with metal to prevent colour development. Warning TMB is toxic avoid direct contact with hands. Dispose of properly
- If a dark blue colour develops within a few minutes after preparation, this indicates that the TMB solution
  has been contaminated and must be discarded. Read absorbance's within 1 hour after completion of the
  assay



- When pipetting reagents, maintain a consistent order of addition from well-to-well. This will ensure equal incubation times for all wells
- · Follow incubation times described in the assay procedure
- · Dispense the TMB solution within 15 min of the washing of the microtitre plate

## 8. Assay Preparation

Bring all reagents to room temperature before use

## 8.1. Assay Design

Determine the number of microwell strips required to test the desired number of samples plus appropriate number of wells needed for running zeros and standards. Each sample, standard and zero should be tested **in duplicate**. Remove sufficient microwell strips for testing from the pouch immediately prior to use. Return any wells not required for this assay with desiccant to the pouch. Seal tightly and return to 2-8°C storage.

**Example plate layout** (example shown for a 6-point standard curve)

	Stand	indards Sample Wells										
	1	2	3	4	5	6	7	8	9	10	11	12
Α	10000	10000										
В	5000	5000										
С	2500	2500										
D	125	125										
Е	625	625										
F	312.5	312.5										
G	Blank	Blank										
Н												

All remaining empty wells can be used to test samples in duplicate

#### 8.2. Preparation of Wash Buffer

Dilute the (200x) wash buffer concentrate 200-fold with distilled water to give a 1x working solution. Pour entire contents (10 ml) of the Washing Buffer Concentrate into a clean 2,000 ml graduated cylinder. Bring final volume to 2,000 ml with glass-distilled or deionized water. Mix gently to avoid foaming. Transfer to a clean wash bottle and store at 2°-8°C for up to 1 week.

#### 8.3. Preparation of Standard Diluent Buffer

Add the contents of the vial (10x concentrate) to 225ml of distilled water before use.

This solution can be stored at 2-8°C for up to 1 week.



#### 8.4. Preparation of Standard

Standard vials must be reconstituted with the volume of standard diluent shown on the vial immediately prior to use. This reconstitution gives a stock solution of 10 000 pg/ml of TRAIL R4. **Mix the reconstituted standard gently by inversion only**. Serial dilutions of the standard are made directly in the assay plate to provide the concentration range from 10 000 to 312.5pg/ml. A fresh standard curve should be produced for each new assay.

- Immediately after reconstitution add 200µl of the reconstituted standard to wells A1 and A2, which provides the highest concentration standard at 10 000 pg/ml
- Add 100µl of Standard Diluent to the remaining standard wells B1 and B2 to F1 and F2
- Transfer 100µl from wells A1 and A2 to B1 and B2. Mix the well contents by repeated aspirations and ejections taking care not to scratch the inner surface of the wells
- Continue this 1:1 dilution using 100µl from wells B1 and B2 through to wells F1 and F2 providing a serial diluted standard curve ranging from 1000pg/ml to 31.25pg/ml
- Discard 100µl from the final wells of the standard curve (F1 and F2)

Alternatively, these dilutions can be performed in separate clean tubes and immediately transferred directly into the relevant wells.

## 8.5. Preparation of Samples

Before testing, human serum or plasmas samples have to be diluted 1:4 in standard buffer diluent.

## 8.6. Preparation of Biotinylated Anti TRAIL R4

It is recommended this reagent is prepared immediately before use. Dilute the biotinylated anti TRAIL R3 with the biotinylated antibody diluent in an appropriate clean glass vial using volumes appropriate to the number of required wells. Please see example volumes below:

Number of wells	Biotinylated	Biotinylated		
required	Antibody (µI)	Antibody Diluent (µI)		
16	40	1060		
24	60	1590		
32	80	2120		
48	120	3180		
96	240	6360		



## 8.7. Preparation of HRP-Conjugate

It is recommended to centrifuge vial for a few seconds in a microcentrifuge to collect all the volume at the bottom.

Dilute the  $5\mu$ l vial with 0.5ml of HRP diluent **immediately before use.** Do-not keep this diluted vial for future experiments. Further dilute the HRP solution to volumes appropriate for the number of required wells in a clean glass vial. Please see example volumes below:

Number of wells required	Streptavidin-HRP (µl)	Streptavidin-HRP Diluent (ml)
16	30	2
24	45	3
32	60	4
48	75	5
96	150	10



## 9. Method

We strongly recommend that every vial is mixed thoroughly without foaming prior to use except the standard vial which must be mixed gently by inversion only.

Prepare all reagents as shown in section 8.

Note: Final preparation of Biotinylated anti TRAIL R4 (section 8.6) and Streptavidin-HRP (section 8.7) should occur immediately before use.

As	ssay Step	Details
1.	Addition	Prepare Standard curve as shown in section 8.4
2.	Addition	Add 100µl of Standard Diluent to <b>zero</b> wells and 100µl of Samples to designated wells
3.	Incubation	Cover with a plastic plate cover and incubate at room temperature (18 to 25°C) for 2 hour(s)
4.	Wash	Remove the cover and wash the plate as follows:  a) Aspirate the liquid from each well b) Dispense 0.3 ml of <b>1x washing solution</b> into each well c) Aspirate the contents of each well d) Repeat step b and c another two times
5.	Addition	Add 50µl of diluted <b>biotinylated anti TRAIL R4</b> to all wells
6.	Incubation	Cover with a plastic plate cover and incubate at room temperature (18 to 25°C) for 1 hour
7.	Wash	Repeat wash step 4.
8.	Addition	Add 100µl of <b>Streptavidin-HRP</b> solution into all wells
9.	Incubation	Cover with a plastic plate cover and incubate at room temperature (18 to 25°C) for <b>30</b> min
10.	Wash	Repeat wash step 4.
11.	Addition	Add 100µl of ready-to-use <b>TMB Substrate Solution</b> into all wells
12.	Incubation	Incubate in the dark for <b>10-20 minutes</b> * at room temperature. Avoid direct exposure to light by wrapping the plate in aluminium foil.
13.	Addition	Add 100µl of <b>H₂SO₄: Stop Reagent</b> into all wells

**Read the absorbance** value of each well (immediately after step 14.) on a spectrophotometer using 450 nm as the primary wavelength and optionally 620 nm as the reference wave length (610 nm to 650 nm is acceptable).

<sup>\*</sup>Incubation time of the substrate solution is usually determined by the ELISA reader performance. Many ELISA readers only record absorbance up to 2.0 O.D. Therefore, the colour development within individual microwells must be observed by the analyst, and the substrate reaction stopped before positive wells are no longer within recordable range



## 10. Data Analysis


Calculate the average absorbance values for each set of duplicate standards and samples. Ideally duplicates should be within 20% of the mean.

Generate a linear standard curve by plotting the average absorbance of each standard on the vertical axis versus the corresponding human TRAIL R4 standard concentration on the horizontal axis.

The amount of TRAIL R4 in each sample is determined by extrapolating OD values against TRAIL R4 standard concentrations using the standard curve.

## **Example TRAIL R4 Standard curve**

Standard	tandard TRAIL R4		CV (%)
1	10 000	1.989	0.6
2	5000	1.142	0.7
3	2500	0.705	7.3
4	1250	0.381	1.3
5	625	0.224	1.7
6	312.5	0.134	9.5
Zero	0	0.039	-

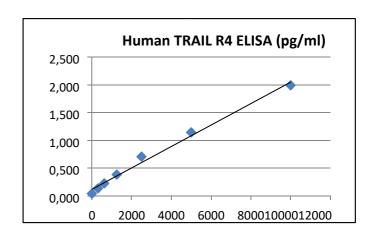

**Note**; curve shown above should not be used to determine results. Every laboratory must produce a standard curve for each set of microwell strips assayed.

For serum and plasma samples which have been diluted according to the protocol (1:4), the calculated concentration should be multiplied by the dilution factor (x4).

# 11. Assay limitations

Do not extrapolate the standard curve beyond the maximum standard curve point. The dose-response is non-linear in this region and good accuracy is difficult to obtain. Concentrated samples above the maximum standard concentration must be diluted with Standard diluent or with your own sample buffer to produce an OD value within the range of the standard curve. Following analysis of such samples always multiply results by the appropriate dilution factor to produce actual final concentration.

The influence of various drugs on end results has not been investigated. Bacterial or fungal contamination and laboratory cross-contamination may also cause irregular results.



Improper or insufficient washing at any stage of the procedure will result in either false positive or false negative results. Completely empty wells before dispensing fresh Washing Buffer, fill with Washing Buffer as indicated for each wash cycle and do not allow wells to sit uncovered or dry for extended periods.

Disposable pipette tips, flasks or glassware are preferred, reusable glassware must be washed and thoroughly rinsed of all detergents before use.

As with most biological assays conditions may vary from assay to assay therefore a fresh standard curve must be prepared and run for every assay.

## 12. Performance Characteristics

## 12.1. Sensitivity

The sensitivity, minimum detectable dose of human TRAIL R4 using this ELISA Genie TRAIL R4 ELISA kit was found to be **64 pg/ml**. This was determined by adding 3 standard deviations to the mean OD obtained when the zero standard was assayed 32 times.

## 12.2. Specificity

The assay recognizes natural and recombinant human TRAIL R4. To define specificity of this ELISA, several proteins were tested for cross reactivity. There was no cross reactivity observed for any protein tested: TRAIL, CD117, IL-6R, IL-2R, CD116, TRAIL R2, TRAIL R3, TRAIL R1, CD178 and Granzyme B.

#### 12.3. Precision

Intra-assay	Inter-assay
-------------	-------------

Sample	n	Mean (pg/mL)	SD	CV%	Sample	n	Mean (pg/mL)	SD	CV%
A	8	5723	212	3.7	A	12	6531	627	9.6
В	8	3063	218	7.1	В	12	3123	369	11.8

#### 12.4. Dilution Parallelism

One human serum containing 10 000 pg/ml of TRAIL R4 was serially diluted in standard buffer diluent. Linear regression of samples versus the expected concentration yielded a correlation coefficient of 0.99.

## 12.5. Spike Recovery

The spike recovery was evaluated by spiking concentrations of natural TRAIL R4 in human serum. Recoveries ranged from 68% to 87% with an overall **mean recovery of 77%.** 

#### 12.6. Expected Serum Values

TRAIL R4 is not detected in healthy donor sera (95 sera tested).



# 13. Assay Summary

Total procedure length: 3h45 mn

Add 100µl of sample or diluted standard

 $\downarrow \! \downarrow$ 

Incubate 2 hours at room temperature

 $\downarrow$ 

Wash three times

 $\downarrow \downarrow$ 

Add 50µl of biotinylated detection antibody

 $\downarrow \downarrow$ 

Incubate 1 hour at room temperature

 $\Downarrow$ 

Wash three times

 $\Downarrow$ 

Add 100µl of Streptavidin-HRP

 $\Downarrow$ 

Incubate 30min at room temperature

 $\Downarrow$ 

Wash three times

 $\Downarrow$ 

Add 100  $\mu I$  of ready-to-use TMB Protect from light. Let the color develop for 10-20 mn.

 $\Downarrow$ 

Add 100 H<sub>2</sub>SO<sub>4</sub>

 $\Downarrow$ 

Read Absorbance at 450 nm



## **Notes**



## <u>Notes</u>



## **Notes**





## **Contact Us:**

Reagent Genie Ireland Limited
The Steelworks
Dublin
Ireland
D01 KA00

Email: info@reagentgenie.com Phone: 00-353-1-887-9802