



Human IL-12p40 + p70 PharmaGenie ELISA Kit

SKU: HUDC0039

Instructions for use

For research use only

Table of Contents

1. Intended use.....	3
2. Introduction.....	3
2.1. Summary.....	3
2.2. Principle of the method.....	3
3. Reagents provided and reconstitution	4
4. Materials required but not provided	4
5. Storage Instructions.....	4
6. Specimen collection, processing & storage	5
7. Safety & precautions for use.....	5
8. Assay Preparation	6
8.1. Assay Design	6
8.2. Preparation of Wash Buffer	6
8.3. Preparation of Standard Diluent Buffer.....	6
8.4. Preparation of Standard	7
8.5. Preparation of Controls	7
8.6. Preparation of Biotinylated anti-IL-12p40	7
8.7. Preparation of Streptavidin-HRP	8
9. Method.....	9
10. Data Analysis	10
11. Assay limitations.....	10
12. Performance Characteristics	10
12.1. Sensitivity.....	10
12.2. Specificity.....	11
12.3. Precision	11
12.4. Dilution Parallelism.....	12
12.5. Spike Recovery	12
12.6. Stability	12
12.7. Expected serum values.....	12
13. Assay Summary	13

Human IL-12p40 + p70 PharmaGenie ELISA Kit

1. Intended use

The ELISA Genie IL-12p40 + p70 PharmaGenie ELISA kit is a solid phase sandwich ELISA for the *in-vitro* qualitative and quantitative determination of the IL-12p40 (p40 subunit) in supernatants, buffered solutions or serum and plasma samples.

This assay will recognise both natural and recombinant human IL-12p40 monomer but also when in homodimer (IL-12p80), heterodimer with the p35 protein (IL-12p70) or with the p19 protein (IL-23).

PharmaGenie ELISA Kits from ELISA Genie are a premium range of pre-coated ELISA kits especially designed for scientists working in pharmaceutical, biotech & CRO sectors. PharmaGenie ELISA kits are produced using high quality monoclonal antibody pairs & optimized reagents that have been manufactured according ISO 9001:2000 quality systems and are excellent assays to help discover our future.

This kit has been configured for research use only. Not suitable for use in therapeutic procedures.

2. Introduction

2.1. Summary

IL-12 is a pleiotropic cytokine initially called cytotoxic lymphocyte maturation factor (CLMF) or natural killer cell stimulatory factor (NKSF) mainly produced by monocytes, macrophages and dendritic cells in response to bacterial products or upon interaction with activated T Cell.

IL-12 induces IFN γ production and increases proliferation and cytotoxic activity of T and NK cells. Moreover, IL-12 induces CD4⁺ polarization to the Th1 phenotype that mediates immunity against intracellular pathogens.

IL-12p70 is the biologically active form and is composed of two subunits: IL-12p35 (also named IL-12A) and IL-12p40 (also named IL-12B). The p35 subunit has homology to IL-6, while p40 has homology with IL-23. The p40 subunit has been found to be expressed in a higher excess over p70. It is present in serum and plasma as a monomer, a heterodimer with IL-12p35 (biologically active IL-12p70) or heterodimer with IL23p19 (IL-23). The subunits are genetically unrelated and are regulated independently: IL-12p40 is produced constitutively and in excess of IL-12p35 and IL-23p19.

Several studies show that IL-12p40, in monomer or homodimer form, was able to bind to IL-12R beta-1 and seem to be antagonist to IL-12p70 and IL-23 binding.

Following these properties and the role of IL-12 in autoimmunity, IL-12p40 was described as a target in the treatment of autoimmune and systemic inflammatory diseases (Crohn's disease, Psoriasis, multiple sclerosis).

2.2. Principle of the method

A capture Antibody highly specific for IL-12p40 has been coated to the wells of the microtiter strip plate provided during manufacture. Binding of IL-12p40 samples and known standards to the capture antibodies and subsequent binding of the biotinylated anti-IL-12p40 secondary antibody to the analyte is completed during the same incubation period. Any excess unbound analyte and secondary antibody is removed. The HRP conjugate solution is then added to every well including the zero wells, following incubation excess conjugate is removed by careful washing.

A chromogen substrate is added to the wells resulting in the progressive development of a blue coloured complex with the conjugate. The colour development is then stopped by the addition of acid turning the

resultant final product yellow. The intensity of the produced coloured complex is directly proportional to the concentration of IL-12p40 present in the samples and standards.

The absorbance of the colour complex is then measured and the generated OD values for each standard are plotted against expected concentration forming a standard curve. This standard curve can then be used to accurately determine the concentration of IL-12p40 in any sample tested.

3. Reagents provided and reconstitution

Reagents (Store @ 2-8°C)	Quantity 1x48 well kit	Quantity 1x96 well kit	Quantity 2x96 well kit	Reconstitution
96 well microtiter strip plate	1/2	1	2	Ready to use (Pre-coated)
Plastic plate covers	2	2	4	n/a
IL-12p40 Standard: 2000 pg/ml	1	2	4	Reconstitute as directed on the vial (see Assay preparation, section 8)
IL-12p40 Control	1	2	4	Reconstitute as directed on the vial (see Assay preparation, section 8)
Standard Diluent (Buffer)	1 (25ml)	1 (25ml)	1 (25ml)	10x Concentrate, dilute in distilled water (see Assay preparation, section 8)
Standard Diluent Serum	1 (7 ml)	1 (7 ml)	2 (7 ml)	Ready to use
Biotinylated antiIL-12p40	1 (0.4ml)	1 (0.4ml)	2 (0.4ml)	Dilute in biotinylated antibody diluent (see Assay preparation, section 8)
Biotinylated Antibody diluent	1 (7ml)	1 (7ml)	1 (13ml)	Ready to use
Streptavidin-HRP	1 (5µl)	2 (5µl)	4 (5µl)	Add 0.5ml of HRP diluent prior to use (see Assay preparation, section 8)
HRP Diluent	1 (23ml)	1 (23ml)	1 (23ml)	Ready to use
Wash Buffer	1 (10ml)	1 (10ml)	2 (10ml)	200x Concentrate dilute in distilled water (see Assay preparation, section 8)
TMB Substrate	1 (11ml)	1 (11ml)	1 (24ml)	Ready to use
H ₂ SO ₄ stop reagent	1 (11ml)	1 (11ml)	2 (11ml)	Ready to use

4. Materials required but not provided

- Microtiter plate reader fitted with appropriate filters (450 nm required with optional 620 nm reference filter)
- Microplate washer or wash bottle
- 10, 50, 100, 200 and 1,000µl adjustable single channel micropipettes with disposable tips
- 50-300µl multi-channel micropipette with disposable tips
- Multichannel micropipette reagent reservoirs
- Distilled water
- Vortex mixer
- Miscellaneous laboratory plastic and/or glass, if possible sterile

5. Storage Instructions

Store kit reagents between 2 and 8°C. Immediately after use remaining reagents should be returned to cold storage (2-8°C). Expiry of the kit and reagents is stated on box front labels. The expiry of the kit components can only be guaranteed if the components are stored properly, and if, in case of repeated use of one component, the reagent is not contaminated by the first handling.

Wash buffer 1X: Once prepared, store at 2-8°C for up to 1 week.

Standard diluent Buffer 1X: Once prepared, store at 2-8°C for up to 1 week.

Reconstituted Standard/Control: Once prepared use immediately and do not store.

Diluted Biotinylated Anti-IL-12p40: Once prepared use immediately and do not store.

Diluted Streptavidin-HRP: Once prepared use immediately and do not store.

6. Specimen collection, processing & storage

Cell culture supernatants, human serum, plasma or other biological samples will be suitable for use in the assay. Remove serum from the clot or red cells, respectively, as soon as possible after clotting and separation.

Cell culture supernatants: Remove particulates and aggregates by spinning at approximately 1000 x g for 10 min.

Serum: Use pyrogen/endotoxin free collecting tubes. Serum should be removed rapidly and carefully from the red cells after clotting. Following clotting, centrifuge at approximately 1000 x g for 10 min and remove serum.

Plasma: EDTA, citrate and heparin plasma can be assayed. Spin samples at 1000 x g for 30 min to remove particulates. Harvest plasma.

Storage: If not analyzed shortly after collection, samples should be aliquoted (250-500µl) to avoid repeated freeze-thaw cycles and stored frozen at -70°C. Avoid multiple freeze-thaw cycles of frozen specimens.

Recommendation: Do not thaw by heating at 37°C or 56°C. Thaw at room temperature and make sure that sample is completely thawed and homogeneous before use. When possible avoid use of badly haemolysed or lipemic sera. If large amounts of particles are present these should be removed prior to use by centrifugation or filtration.

7. Safety & precautions for use

- Handling of reagents, serum or plasma specimens should be in accordance with local safety procedures, e.g. CDC/NIH Health manual: "Biosafety in Microbiological and Biomedical Laboratories" 1984. Laboratory gloves should be worn at all times.
- Avoid any skin contact with H₂SO₄ and TMB. In case of contact, wash thoroughly with water.
- Do not eat, drink, smoke or apply cosmetics where kit reagents are used.
- Do not pipette by mouth.
- When not in use, kit components should be stored refrigerated as indicated on vials or bottles labels.
- All reagents should be warmed to room temperature before use. Lyophilized standards should be discarded after use.
- Once the desired number of strips has been removed, immediately reseal the bag to protect the remaining strips from deterioration.
- Cover or cap all reagents when not in use.
- Do not mix or interchange reagents between different lots.
- Do not use reagents beyond the expiration date of the kit.
- Use a clean disposable plastic pipette tip for each reagent, standard, or specimen addition in order to avoid cross contamination, for the dispensing of H₂SO₄ and substrate solution, avoid pipettes with metal parts.
- Use a clean plastic container to prepare the washing solution.
- Thoroughly mix the reagents and samples before use by agitation or swirling.
- All residual washing liquid must be drained from the wells by efficient aspiration or by decantation followed by tapping the plate forcefully on absorbent paper. Never insert absorbent paper directly into the wells.
- The TMB solution is light sensitive. Avoid prolonged exposure to light. Also, avoid contact of the TMB solution with metal to prevent colour development. Warning TMB is toxic avoid direct contact with hands. Dispose of properly.

- If a dark blue colour develops within a few minutes after preparation, this indicates that the TMB solution has been contaminated and must be discarded. Read absorbance's within 1 hour after completion of the assay.
- When pipetting reagents, maintain a consistent order of addition from well-to-well. This will ensure equal incubation times for all wells.
- Follow incubation times described in the assay procedure.
- Dispense the TMB solution within 15 min of the washing of the microtiter plate.

8. Assay Preparation

Bring all reagents to room temperature before use

8.1. Assay Design

Determine the number of microwell strips required to test the desired number of samples plus appropriate number of wells needed for running zeros and standards. Each sample, standard, zero and control should be tested **in duplicate**. Remove sufficient microwell strips for testing from the pouch immediately prior to use. Return any wells not required for this assay with desiccant to the pouch. Seal tightly and return to 28°C storage.

Example plate layout (example shown for a 6-point standard curve)

	Standards / Controls		Sample Wells									
	1	2	3	4	5	6	7	8	9	10	11	12
A	2000	2000										
B	1000	1000										
C	500	500										
D	250	250										
E	125	125										
F	62.5	62.5										
G	zero	zero										
H	Ctrl	Ctrl										

All remaining empty wells can be used to test samples in duplicate

8.2. Preparation of Wash Buffer

If crystals have formed in the concentrate Wash buffer, warm it gently until complete dissolution. Dilute the (200X) wash buffer concentrate 200-fold with distilled water to give a 1X working solution. Pour entire contents (10 ml) of the Concentrate Washing Buffer into a clean 2,000 ml graduated cylinder. Bring final volume to 2,000 ml with glass-distilled or deionized water. Mix gently to avoid foaming. Transfer to a clean wash bottle and store at 2°-25°C.

8.3. Preparation of Standard Diluent Buffer

If crystals have formed in the concentrate Standard Diluent Buffer, warm it gently until complete dissolution. Dilute the (10X) standard diluent concentrate 10-fold with distilled water to give a 1X working solution. Pour entire contents (25 ml) of the Concentrate Standard diluent into a clean 250 ml graduated cylinder. Bring final volume to 250 ml with glass-distilled or deionized water. Transfer to a clean wash bottle and store at 2°-25°C.

This solution can be stored at 2-8°C for up to 1 week.

8.4. Preparation of Standard

Depending on the type of samples you are assaying, the kit may include two standard diluents. Because biological fluids might contain proteases or cytokine-binding proteins that could modify the recognition of the cytokine you want to measure, you should reconstitute standard vials with the most appropriate Standard Diluent.

For **serum and plasma** samples: use Standard Diluent - Serum For
cell culture supernatants: use Standard Diluent Buffer

Standard vials must be reconstituted with the volume of standard diluent shown on the vial immediately prior to use. This reconstitution gives a stock solution of 2000 pg/ml of IL-12p40. Mix the reconstituted standard gently by inversion only. Serial dilutions of the standard are made directly in the assay plate to provide the concentration range from 2000 to 62.5 pg/ml. A fresh standard curve should be produced for each new assay.

- Immediately after reconstitution add 200µl of the reconstituted standard to wells A1 and A2, which provides the highest concentration standard at 2000 pg/ml.
- Add 100µl of standard diluent to the remaining standard wells B1 and B2 to F1 and F2.
- Transfer 100µl from wells A1 and A2 to B1 and B2. Mix the well contents by repeated aspirations and ejections taking care not to scratch the inner surface of the wells.
- Continue this 1:1 dilution using 100µl from wells B1 and B2 through to wells F1 and F2 providing a serial diluted standard curve ranging from 2000 pg/ml to 62.5 pg/ml.
- Discard 100µl from the final wells of the standard curve (F1 and F2).

Alternatively, these dilutions can be performed in separate clean tubes and immediately transferred into the relevant wells.

8.5. Preparation of Controls

Freeze-dried control vials should also be reconstituted with the most appropriate Standard Diluent to your samples.

For **serum and plasma** samples: use Standard Diluent - Serum For
cells culture supernatants: use Standard Diluent Buffer

The supplied Controls must be reconstituted with the volume of Standard Diluent indicated on the vial. Reconstitution of the freeze-dried material with the indicated volume, will give a solution at the concentration stated on the vial. Do not store after use.

8.6. Preparation of Biotinylated anti-IL-12p40

It is recommended this reagent is prepared immediately before use. Dilute the biotinylated anti-IL-12p40 with the biotinylated antibody diluent in an appropriate clean glass vial using volumes appropriate to the number of required wells. Please see example volumes below:

Number of wells required	Biotinylated Antibody (µl)	Biotinylated Antibody Diluent (µl)
16	40	1060
24	60	1590
32	80	2120
48	120	3180
96	240	6360

8.7. Preparation of Streptavidin-HRP

It is recommended to centrifuge vial for a few seconds in a microcentrifuge to collect all the volume at the bottom.

Dilute the 5µl vial with 0.5ml of HRP diluent **immediately before use**. Do not keep this diluted vial for future experiments. Further dilute the HRP solution to volumes appropriate for the number of required wells in a clean glass vial. Please see example volumes below:

Number of wells required	Streptavidin-HRP (µl)	Streptavidin-HRP Diluent (ml)
16	30	2
24	45	3
32	60	4
48	75	5
96	150	10

9. Method

We strongly recommend that every vial is mixed thoroughly without foaming prior to use.

Prepare all reagents as shown in section 8.

Note: final preparation of Biotinylated Secondary Antibody (section 8.6) and Streptavidin-HRP (section 8.7) should occur immediately before use.

Assay Step		Details
1.	Addition	Prepare Standard curve as shown in section 8.4 above and add in duplicate to appropriate wells
2.	Addition	Add 100µl of each Sample, Control and zero (Standard diluent) in duplicate to appropriate number of wells
3.	Addition	Add 50µl of diluted biotinylated anti-IL-12p40 to all wells
4.	Incubation	Cover with a plastic plate cover and incubate at room temperature (18 to 25°C) for 1 hour
5.	Wash	Remove the cover and wash the plate as follows: a) Aspirate the liquid from each well b) Dispense 0.3 ml of 1x washing solution into each well c) Aspirate the contents of each well d) Repeat step b and c another two times
6.	Addition	Add 100µl of Streptavidin-HRP solution into all wells
7.	Incubation	Cover with a plastic plate cover and incubate at room temperature (18 to 25°C) for 30 min
8.	Wash	Repeat wash step 5.
9.	Addition	Add 100µl of ready-to-use TMB Substrate Solution into all wells
10.	Incubation	Incubate in the dark for 12-15 minutes* at room temperature. Avoid direct exposure to light by wrapping the plate in aluminium foil.
11.	Addition	Add 100µl of H₂SO₄: Stop Reagent into all wells
Read the absorbance value of each well (immediately after step 11.) on a spectrophotometer using 450 nm as the primary wavelength and optionally 620 nm as the reference wave length (610 nm to 650 nm is acceptable).		

** Incubation time of the substrate solution is usually determined by the ELISA reader performance. Many ELISA readers only record absorbance up to 2.0 O.D. Therefore, the colour development within individual microwells must be observed by the analyst, and the substrate reaction stopped before positive wells are no longer within recordable range.*

10. Data Analysis


Calculate the average absorbance values for each set of duplicate standards, controls and samples. Ideally duplicates should be within 20% of the mean.

Generate a linear standard curve by plotting the average absorbance of each standard on the vertical axis versus the corresponding IL-12p40 standard concentration on the horizontal axis.

The amount of IL-12p40 in each sample is determined by extrapolating OD values against IL-12p40 standard concentrations using the standard curve.

Example IL-12p40 Standard curve

Standard	IL-12p40 Concentration (pg/ml)	OD (450nm) mean	CV (%)
1	2000	2.968	1.64
2	1000	1.737	7.16
3	500	0.942	6.15
4	250	0.536	1.65
5	125	0.317	5.95
6	62.5	0.196	12.28
zero	0	0.063	38.81

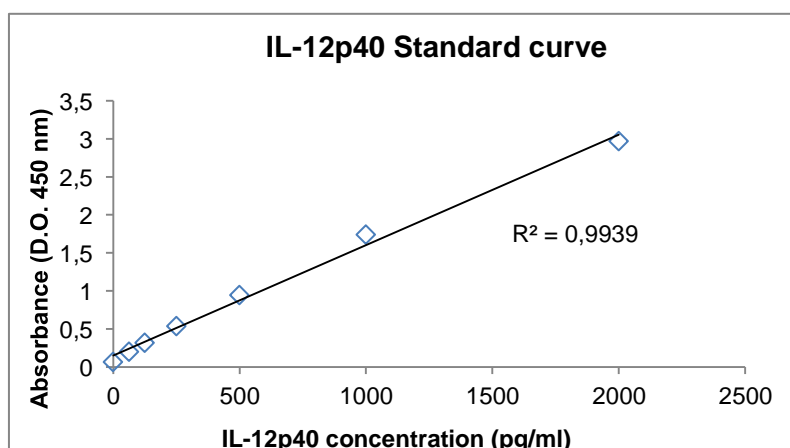

Note: curve shown above should not be used to determine results. Every laboratory must produce a standard curve for each set of microwell strips assayed.

11. Assay limitations

Do not extrapolate the standard curve beyond the maximum standard curve point. The dose-response is non-linear in this region and good accuracy is difficult to obtain. Concentrated samples above the maximum standard concentration must be diluted with Standard diluent or with your own sample buffer to produce an OD value within the range of the standard curve. Following analysis of such samples always multiply results by the appropriate dilution factor to produce actual final concentration.

The influence of various drugs on end results has not been investigated. Bacterial or fungal contamination and laboratory cross-contamination may also cause irregular results.

Improper or insufficient washing at any stage of the procedure will result in either false positive or false negative results. Completely empty wells before dispensing fresh Wash Buffer, fill with Wash Buffer as indicated for each wash cycle and do not allow wells to sit uncovered or dry for extended periods.

Disposable pipette tips, flasks or glassware are preferred, reusable glassware must be washed and thoroughly rinsed of all detergents before use.

As with most biological assays conditions may vary from assay to assay therefore **a fresh standard curve must be prepared and run for every assay.**

12. Performance Characteristics

12.1. Sensitivity

The sensitivity or minimum detectable dose of IL-12p40 using this ELISA Genie IL-12p40 ELISA kit was found to be **20 pg/ml**. This was determined by adding 3 standard deviations to the mean OD obtained when the zero standard was assayed 40 times.

12.2. Specificity

The assay recognizes both natural and recombinant human IL-12p40. To define the specificity of this ELISA several proteins were tested for cross reactivity. There was no cross reactivity observed for any protein tested: IL-1 α , IL-1 β , IL-10, IL-2, IFN γ , IL-4, IL-6, TNF α , IL-8 and IL-13.

The assay will also recognize the p40 when in heterodimer with the p35 protein (IL-12p70) or with the p19 protein (IL-23).

12.3. Precision

Intra-assay

Reproducibility within the assay was evaluated in three independent experiments. Each assay was carried out with 6 replicates (3 duplicates) of samples containing different concentrations of IL-12p40: 2 in human pooled serum and 3 in standard diluent. 1 standard curve was run on each plate. Data below show the mean IL-12p40 concentration and the coefficient of variation for each sample. **The calculated overall coefficient of variation was 3.7%.**

Session	Sample	Mean IL-12p40 pg/ml	SD	CV%
Session 1	Sample 1	2010.33	67.86	3.4
	Sample 2	996.67	32.32	3.2
	Sample 3	474.00	4.58	1.0
	Sample 4	1314.00	59.81	4.6
	Sample 5	687.67	31.18	4.5
Session 2	Sample 1	2345.00	33.29	1.4
	Sample 2	1254.33	62.98	5.0
	Sample 3	511.00	11.36	2.2
	Sample 4	1579.00	61.29	3.9
	Sample 5	842.67	49.50	5.9
Session 3	Sample 1	2206.00	80.58	3.7
	Sample 2	1155.00	22.54	2.0
	Sample 3	508.00	29.82	5.9
	Sample 4	1456.67	55.37	3.8
	Sample 5	825.00	47.51	5.8

Inter-assay

Assay to assay reproducibility within one laboratory was evaluated in three independent experiments by two technicians. Each assay was carried out with 6 replicates (3 duplicates) of samples containing different concentrations of IL-12p40: 2 in human pooled serum and 3 in standard diluent. 1 standard curve was run on each plate. Data below show the mean IL-12p40 concentration and the coefficient of variation for each sample.

The calculated overall coefficient of variation was 9.5%.

	Sample 1	Sample 2	Sample 3	Sample 4	Sample 5
Mean IL-12p40 pg/ml	2277	1199	479	1488	826
SD	160	108	49	123	106
CV%	7.0	9.0	10.3	8.2	12.8

12.4. Dilution Parallelism

In two independent experiments two spiked human serum samples with different levels of IL-12p40 were analysed at different serial two-fold dilutions (1:2 to 1:8) with two replicates each. Recoveries ranged from 84 to 118% with an overall **mean recovery of 98%**.

12.5. Spike Recovery

The spike recovery was evaluated by spiking 3 concentrations of IL-12p40 in human serum and culture medium in 3 separate experiments.

Recoveries ranged from 99 to 110% with an overall **mean recovery of 103%**.

12.6. Stability

Storage Stability

Aliquots of spiked serum and spiked medium were stored at -20°C , $+2-8^{\circ}\text{C}$, room temperature (RT) and at 37°C and the IL-12p40 level determined after 24h.

There was no significant loss of IL-12p40 reactivity in medium during storage at $+2-8^{\circ}\text{C}$, RT and 37°C . However, a significant loss was observed in serum during storage at $+2-8^{\circ}\text{C}$ (18%), RT (29%) and 37°C (24%).

Freeze-thaw Stability

Aliquots of spiked serum and spiked medium were stored frozen at -20°C and thawed up to 5 times and the IL-12p40 level was determined. There was a slight loss of IL-12p40 after 5 cycles of freezing and thawing.

12.7. Expected serum values

A panel of 24 human sera samples were tested for IL-12p40. See results below:

Number of samples evaluated	Range (pg/ml)	Mean (pg/ml)	Standard deviation (pg/ml)
3	Below the detection level (20pg/ml)		
21	[65;483]	187	126

13. Assay Summary

Total procedure length: 1h45min

**Add 100µl of sample, controls and diluted standard
and 50µl Biotinylated anti-IL-12p40**



Incubate 1 hour at room temperature



Wash three times



Add 100µl of Streptavidin-HRP



Incubate 30 min at room temperature



Wash three times



Add 100µl of ready-to-use TMB
Protect from light. Let the color develop for 12-15 min.



Add 100µl H₂SO₄



Read Absorbance at 450 nm

Notes

Notes



Contact Us:

Reagent Genie Ireland Limited
The Steelworks
Dublin
Ireland
D01 KA00

Email: info@reagentgenie.com

Phone: 00-353-1-887-9802