

## MMP2 Antibody

CAB11144

### Description

---

This MMP2 Antibody is supplied as a kit for advanced applications. The kit includes Bradford Reagent to quantify total protein concentration for accurate sample normalization (Optional).

### Product Information

---

|                      |                                                           |
|----------------------|-----------------------------------------------------------|
| <b>SKU:</b>          | CAB11144                                                  |
| <b>Contents:</b>     | 20 $\mu$ L, 100 $\mu$ L<br>Bradford Reagent: 1 vial (2ml) |
| <b>Category:</b>     | Polyclonal Antibody                                       |
| <b>Synonyms:</b>     | CLG4, MONA, CLG4A, MMP-2, TBE-1, MMP-II, MMP2             |
| <b>Clone:</b>        | -                                                         |
| <b>Applications:</b> | <span>WB</span> <span>IF/ICC</span> <span>ELISA</span>    |
| <b>Conjugation:</b>  | Unconjugated                                              |
| <b>Reactivity:</b>   | Human, Mouse, Rat                                         |

### Antibody Data

---

|                       |                       |
|-----------------------|-----------------------|
| <b>Gene ID:</b>       | 4313                  |
| <b>Uniprot:</b>       | AB_2758429            |
| <b>Host Species:</b>  | Rabbit                |
| <b>Purification:</b>  | Affinity purification |
| <b>Observed MW:</b>   | 72kDa                 |
| <b>Calculated MW:</b> | 74kDa                 |

## Preparation & Storage

**Storage:** Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS containing 50% glycerol, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.

Store Bradford Reagent at Room Temperature for 1 Year.


**Positive Sample:** Mouse lung. Rat lung

**Recommended Dilutions:**


|               |                                                                                                                             |
|---------------|-----------------------------------------------------------------------------------------------------------------------------|
| <b>WB</b>     | 1:500 - 1:1000                                                                                                              |
| <b>IF/ICC</b> | 1:50 - 1:200                                                                                                                |
| <b>ELISA</b>  | Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements. |

**Protein Quantification (Optional):** To quantify total protein levels, use the Bradford Reagent included in this kit. Visit <https://www.assaygenie.com/bradford-protein-assay-protocol/> to view the full protocol

## Validation Data

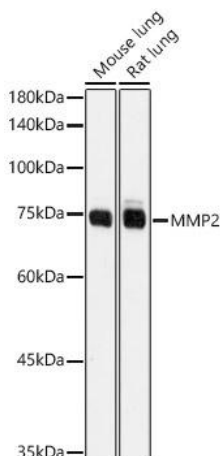

Western blot analysis of various lysates, using Rabbit pAb (CAB11144) at 1:1000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (CABS014) at 1:10000 dilution. Lysates/proteins: 25µg per lane. Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (AbGn00020). Exposure time: 1s.

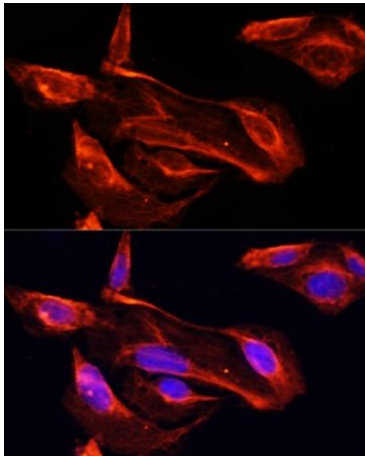

Immunofluorescence analysis of U2OS cells using Rabbit pAb (CAB11144) at dilution of 1:100. Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (CABS007) at 1:500 dilution. Blue: DAPI for nuclear staining.