

## Asymmetric DiMethyl-Histone H3-R8 Antibody

CAB3157

### Description

---

This Asymmetric DiMethyl-Histone H3-R8 Antibody is supplied as a kit for advanced applications. The kit includes Bradford Reagent to quantify total protein concentration for accurate sample normalization (Optional).

### Product Information

---

<b>SKU:</b>	CAB3157
<b>Contents:</b>	20 µL, 100 µL Bradford Reagent: 1 vial (2ml)
<b>Category:</b>	Polyclonal Antibody
<b>Synonyms:</b>	H3t, H3.4, H3/g, H3FT, H3C16, HIST3H3, Asymmetric DiMethyl-Histone H3-R8
<b>Clone:</b>	-
<b>Applications:</b>	WB IHC-P IF/ICC ChIP ELISA
<b>Conjugation:</b>	Unconjugated
<b>Reactivity:</b>	Human, Mouse, Rat, Other (Wide Range Predicted)

### Antibody Data

---

<b>Gene ID:</b>	8290 8350
<b>Uniprot:</b>	AB_2764951
<b>Host Species:</b>	Rabbit
<b>Purification:</b>	Affinity purification
<b>Observed MW:</b>	17 kDa
<b>Calculated MW:</b>	15 kDa

## Preparation & Storage

**Storage:** Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS containing 50% glycerol, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.

Store Bradford Reagent at Room Temperature for 1 Year.


**Positive Sample:** HeLa, NIH/3T3

**Recommended Dilutions:**


<b>WB</b>	1:100 - 1:500
<b>IHC-P</b>	1:50 - 1:200
<b>IF/ICC</b>	1:50 - 1:200
<b>ELISA</b>	Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements. ChIP 5µg antibody for 5µg-10µg of Chromatin

**Protein Quantification (Optional):** To quantify total protein levels, use the Bradford Reagent included in this kit. Visit <https://www.assaygenie.com/bradford-protein-assay-protocol/> to view the full protocol


## Validation Data

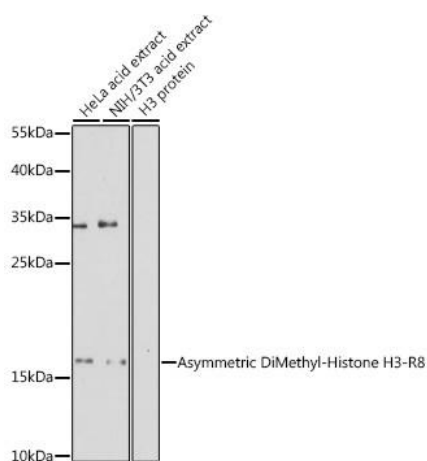


Western blot analysis of various lysates using Asymmetric DiMethyl-Histone - Rabbit pAb (CAB3157) at 1:500 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (CABS014) at 1:10000 dilution. Lysates/proteins: 25µg per lane. Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (AbGn00020). Exposure time: 180s.

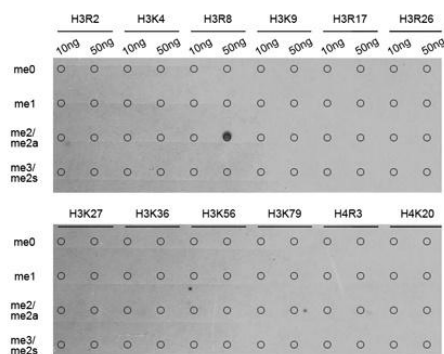


Dot-blot analysis of all sorts of methylation peptides using Asymmetric DiMethyl-Histone - antibody (CAB3157).

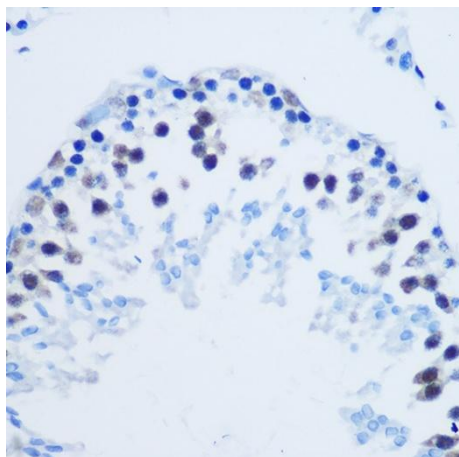

Immunohistochemistry analysis of paraffin-embedded Rat testis using Asymmetric DiMethyl-Histone - Rabbit pAb (CAB3157) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

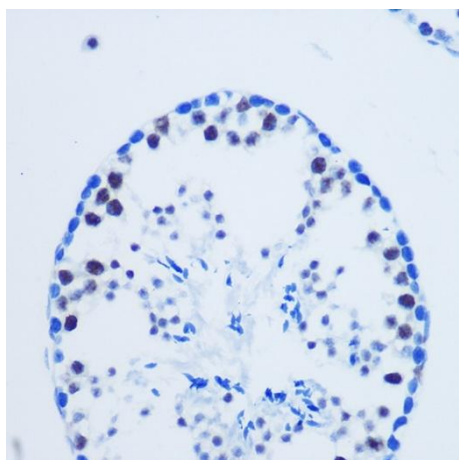

Immunohistochemistry analysis of paraffin-embedded Mouse testis using Asymmetric DiMethyl-Histone - Rabbit pAb (CAB3157) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

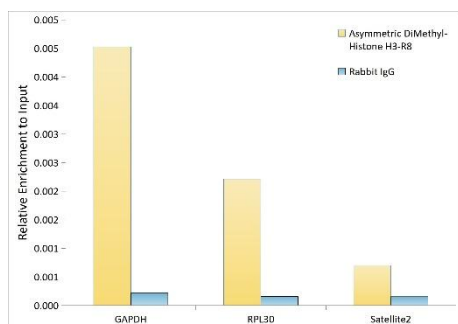

Chromatin immunoprecipitation analysis of extracts of HeLa cells, using Asymmetric DiMethyl-Histone - antibody (CAB3157) and rabbit IgG. The amount of immunoprecipitated DNA was checked by quantitative PCR. Histogram was constructed by the ratios of the immunoprecipitated DNA to the input.