

## MYBPC2 Antibody

CAB6934

### Description

---

This MYBPC2 Antibody is supplied as a kit for advanced applications. The kit includes Bradford Reagent to quantify total protein concentration for accurate sample normalization (Optional).

### Product Information

---

|                      |                                                       |
|----------------------|-------------------------------------------------------|
| <b>SKU:</b>          | CAB6934                                               |
| <b>Contents:</b>     | 20 µL, 100 µL<br>Bradford Reagent: 1 vial (2ml)       |
| <b>Category:</b>     | Polyclonal Antibody                                   |
| <b>Synonyms:</b>     | MYBPC, MYBPCF, fsMyBP-C, MYBPC2                       |
| <b>Clone:</b>        | -                                                     |
| <b>Applications:</b> | <span>WB</span> <span>IHC-P</span> <span>ELISA</span> |
| <b>Conjugation:</b>  | Unconjugated                                          |
| <b>Reactivity:</b>   | Mouse, Rat                                            |

### Antibody Data

---

|                       |                       |
|-----------------------|-----------------------|
| <b>Gene ID:</b>       | 4606                  |
| <b>Uniprot:</b>       | AB_2767492            |
| <b>Host Species:</b>  | Rabbit                |
| <b>Purification:</b>  | Affinity purification |
| <b>Observed MW:</b>   | 135kDa                |
| <b>Calculated MW:</b> | 128kDa                |

## Preparation & Storage

**Storage:** Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS containing 50% glycerol, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.


Store Bradford Reagent at Room Temperature for 1 Year.

**Positive Sample:** Mouse skeletal muscle, Rat skeletal muscle

|                               |              |                                                                                                                             |
|-------------------------------|--------------|-----------------------------------------------------------------------------------------------------------------------------|
| <b>Recommended Dilutions:</b> | <b>WB</b>    | 1:500 - 1:2000                                                                                                              |
|                               | <b>IHC-P</b> | 1:50 - 1:200                                                                                                                |
|                               | <b>ELISA</b> | Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements. |

**Protein Quantification (Optional):** To quantify total protein levels, use the Bradford Reagent included in this kit. Visit <https://www.assaygenie.com/bradford-protein-assay-protocol/> to view the full protocol

## Validation Data

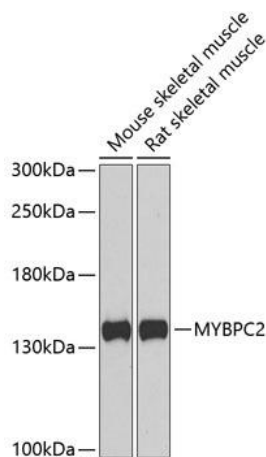

Western blot analysis of various lysates using MYBPC2 Rabbit pAb (CAB6934) at 1:1000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (CABS014) at 1:10000 dilution. Lysates/proteins: 25µg per lane. Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (AbGn00020). Exposure time: 90s.