

Proprotein Convertase 9(PCSK9) Antibody

CAB7860

Description

This Proprotein Convertase 9(PCSK9) Antibody is supplied as a kit for advanced applications. The kit includes Bradford Reagent to quantify total protein concentration for accurate sample normalization (Optional).

Product Information

SKU:	CAB7860
Contents:	20 µL, 100 µL Bradford Reagent: 1 vial (2ml)
Category:	Polyclonal Antibody
Synonyms:	FH3, PC9, FHCL3, NARC1, LDLCQ1, NARC-1, HCHOLA3, Proprotein Convertase 9(PCSK9)
Clone:	-
Applications:	WB IHC-P IF/ICC IP ELISA
Conjugation:	Unconjugated
Reactivity:	Human, Mouse, Rat

Antibody Data

Gene ID:	255738
Uniprot:	AB_2770818
Host Species:	Rabbit
Purification:	Affinity purification
Observed MW:	80 kDa/65 kDa
Calculated MW:	74 kDa

Preparation & Storage

Storage: Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS containing 50% glycerol, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.

Store Bradford Reagent at Room Temperature for 1 Year.


Positive Sample: HeLa, Mouse liver, Rat liver, Rat testis

Recommended Dilutions:


WB	1:500 - 1:5000
IP	0.5µg-4µg antibody for 200µg-400µg extracts of whole cells
IF/ICC	1:50 - 1:200
IHC-P	1:50 - 1:200
ELISA	Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

Protein Quantification (Optional): To quantify total protein levels, use the Bradford Reagent included in this kit. Visit <https://www.assaygenie.com/bradford-protein-assay-protocol/> to view the full protocol

Validation Data

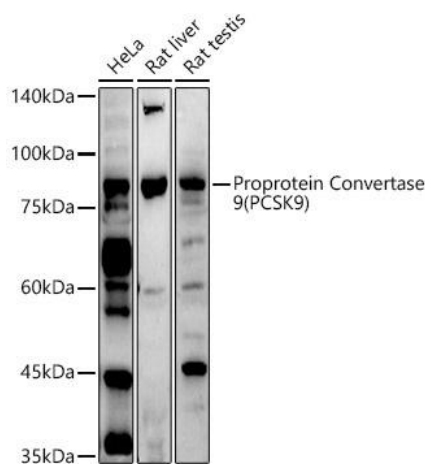


Western blot analysis of various lysates using Proprotein Convertase 9(PCSK9) Rabbit pAb (CAB7860) at 1:1000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (CABS014) at 1:10000 dilution. Lysates/proteins: 25µg per lane. Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (AbGn00020). Exposure time: 180s.

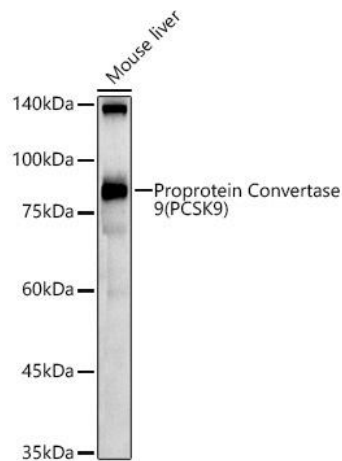


Western blot analysis of lysates from Mouse liver, using Proprotein Convertase 9(PCSK9) Rabbit pAb (CAB7860) at 1:1000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (CABS014) at 1:10000 dilution. Lysates/proteins: 25µg per lane. Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Enhanced Kit (AbGn00021). Exposure time: 30s.




Immunohistochemistry analysis of paraffin-embedded Mouse intestine tissue using Proprotein Convertase 9(PCSK9) Rabbit pAb (CAB7860) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

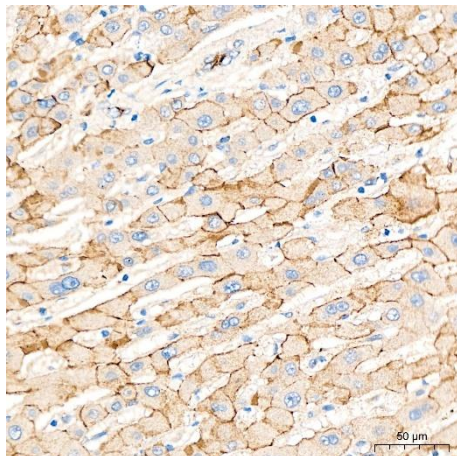


Immunohistochemistry analysis of paraffin-embedded Human liver cancer tissue using Proprotein Convertase 9(PCSK9) Rabbit pAb (CAB7860) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

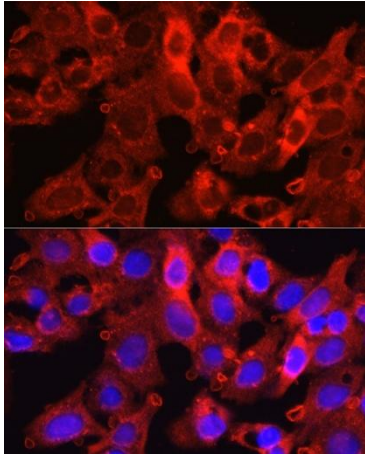


Immunofluorescence analysis of HepG2 cells using Proprotein Convertase 9(PCSK9) Rabbit pAb (CAB7860) at dilution of 1:200 (40x lens). Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (CABS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

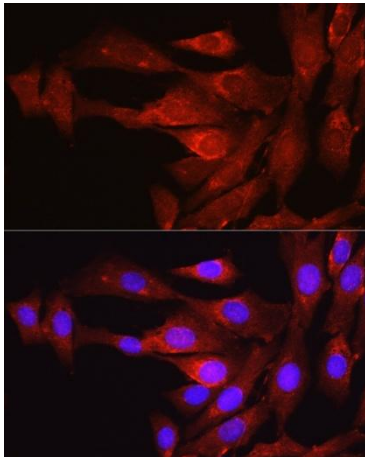

Immunofluorescence analysis of NIH/3T3 cells using Proprotein Convertase 9(PCSK9) Rabbit pAb (CAB7860) at dilution of 1:200 (40x lens). Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (CABS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

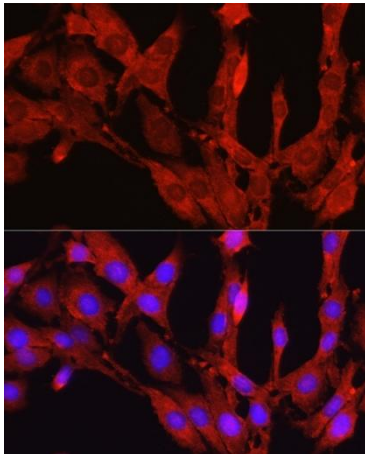

Immunofluorescence analysis of PC-12 cells using Proprotein Convertase 9(PCSK9) Rabbit pAb (CAB7860) at dilution of 1:200 (40x lens). Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (CABS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

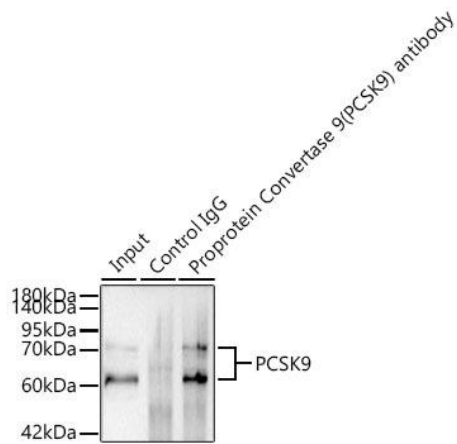

Immunoprecipitation of Proprotein Convertase 9(PCSK9) from 200 µg extracts of HeLa cells was performed using 0.5 µg of Proprotein Convertase 9(PCSK9) Rabbit pAb (CAB7860). Rabbit IgG isotype control (CABC042) was used to precipitate the Control IgG sample. IP samples were eluted with 1X Laemmli Buffer. The Input lane represents 10% of the total input. Western blot analysis of immunoprecipitates was conducted using Proprotein Convertase 9(PCSK9) Rabbit pAb (CAB7860) at a dilution of 1:1000.