

HTR5B Antibody

PACO02280

Description

This HTR5B Antibody is supplied as a kit for advanced applications. The kit includes Bradford Reagent to quantify total protein concentration for accurate sample normalization (Optional).

Product Information

SKU:	PACO02280
Contents:	50µg Bradford Reagent: 1 vial (2ml)
Category:	-
Synonyms:	Htr5b, 5ht5b, 5-hydroxytryptamine receptor 5B, 5-HT-5B, 5-HT5B, Serotonin receptor 5B
Clone:	-
Applications:	WB IF ELISA
Conjugation:	Non-conjugated
Reactivity:	Mouse,Rat

Antibody Data

Isotype:	IgG
Uniprot:	P31387
Host Species:	Rabbit
Purification:	The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen.
Immunogen:	Synthesized peptide derived from the C-terminal region of Human SR-5B.
Immunogen Species:	Mus musculus (Mouse)
Buffer:	Liquid in PBS containing 50% glycerol, 0.5% BSA and 0.02% sodium azide.
Form:	Liquid

Preparation & Storage


Storage: Upon receipt, store at -20°C or -80°C. Avoid repeated freeze.
Store Bradford Reagent at Room Temperature for 1 Year.

Recommended Dilutions:	Application	Recommended Dilution
	WB	1:500-1:2000
	IF	1:200-1:1000
	ELISA	1:40000

Protein Quantification (Optional): To quantify total protein levels, use the Bradford Reagent included in this kit. Visit <https://www.assaygenie.com/bradford-protein-assay-protocol/> to view the full protocol

Validation Data

Image

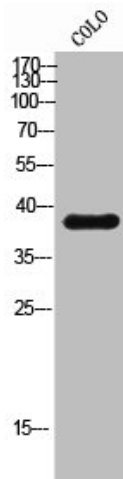

Description

Western Blot analysis of COLO cells using SR-5B Polyclonal Antibody