

KLHL3 Antibody

PACO23772

Description

This KLHL3 Antibody is supplied as a kit for advanced applications. The kit includes Bradford Reagent to quantify total protein concentration for accurate sample normalization (Optional).

Product Information

SKU:	PACO23772
Contents:	100µl Bradford Reagent: 1 vial (2ml)
Category:	Polyclonal Antibody
Synonyms:	KLHL3 antibody, KIAA1129 antibody, Kelch-like protein 3 antibody
Clone:	Polyclonal
Applications:	ELISA WB IHC IF
Conjugation:	Non-conjugated
Reactivity:	Human

Antibody Data

Isotype:	-
Uniprot:	Q9UH77
Host Species:	Rabbit
Purification:	The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen.
Immunogen:	Synthesized peptide derived from N-terminal of Human KLHL3.
Immunogen Species:	Homo sapiens (Human)
Buffer:	-
Form:	Rabbit IgG in phosphate buffered saline (without Mg ²⁺ and Ca ²⁺), pH 7.4, 150mM NaCl, 0.02% sodium azide and 50% glycerol.

Manufacturers Statement: This final kit system is assembled and quality-released by Assay Genie Limited.

Preparation & Storage

Storage: Upon receipt, store at -20°C or -80°C. Avoid repeated freeze.
Store Bradford Reagent at Room Temperature for 1 Year.


Recommended Dilutions:

Application	Recommended Dilution
WB	1:500-1:3000
IHC	1:50-1:100
IF	1:100-1:500

Protein Quantification (Optional): To quantify total protein levels, use the Bradford Reagent included in this kit. Visit <https://www.assaygenie.com/bradford-protein-assay-protocol/> to view the full protocol

Validation Data

Image

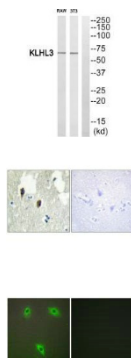

Description

Western blot analysis of extracts from RAW264.7 cells and 3T3 cells, using KLHL3 antibody.

Immunohistochemistry analysis of paraffin-embedded human brain tissue, using KLHL3 antibody.

Immunofluorescence analysis of A549 cells, using KLHL3 antibody.