

ndr1 Antibody

PACO61770

Description

This ndr1 Antibody is supplied as a kit for advanced applications. The kit includes Bradford Reagent to quantify total protein concentration for accurate sample normalization (Optional).

Product Information

| | |
|----------------------|--|
| SKU: | PACO61770 |
| Contents: | 50µl Bradford Reagent: 1 vial (2ml) |
| Category: | - |
| Synonyms: | - |
| Clone: | Polyclonal |
| Applications: | ELISA WB |
| Conjugation: | Non-conjugated |
| Reactivity: | Zebrafish |

Antibody Data

| | |
|---------------------------|---|
| Isotype: | IgG |
| Uniprot: | F1QDT8 |
| Host Species: | Rabbit |
| Purification: | >95%, Protein G purified |
| Immunogen: | Recombinant Zebrafish Nodal-related 1 protein (71-213AA) |
| Immunogen Species: | Zebrafish |
| Buffer: | Preservative: 0.03% Proclin 300 Constituents: 50% Glycerol, 0.01M PBS, pH 7.4 |
| Form: | Liquid |

Preparation & Storage


Storage: Upon receipt, store at -20°C or -80°C. Avoid repeated freeze.
Store Bradford Reagent at Room Temperature for 1 Year.

| Recommended Dilutions: | Application | Recommended Dilution |
|------------------------|-------------|----------------------|
| | WB | 1:500-1:5000 |

Protein Quantification (Optional): To quantify total protein levels, use the Bradford Reagent included in this kit. Visit <https://www.assaygenie.com/bradford-protein-assay-protocol/> to view the full protocol

Validation Data

Image

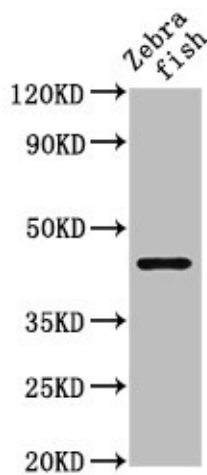

Description

Western Blot Positive WB detected in: Zebrafish tissue
All lanes: ndr1 antibody at 5µg/ml Secondary Goat polyclonal to rabbit IgG at 1/50000 dilution Predicted band size: 45 kDa Observed band size: 45 kDa