

Nampt Antibody

PACO64323

Description

This Nampt Antibody is supplied as a kit for advanced applications. The kit includes Bradford Reagent to quantify total protein concentration for accurate sample normalization (Optional).

Product Information

SKU:	PACO64323
Contents:	50µl Bradford Reagent: 1 vial (2ml)
Category:	-
Synonyms:	Nampt antibody, Pbef1 antibody, Nicotinamide phosphoribosyltransferase antibody, NAMPRTase antibody, Nampt antibody, EC 2.4.2.12 antibody, Pre-B-cell colony-enhancing factor 1 homolog antibody, PBEF antibody, Visfatin antibody
Clone:	Polyclonal
Applications:	ELISA WB
Conjugation:	Non-conjugated
Reactivity:	Rat, Human

Antibody Data

Isotype:	IgG
Uniprot:	Q80Z29
Host Species:	Rabbit
Purification:	>95%, Antigen Affinity purified & Affinity purified
Immunogen:	Recombinant Rat Nampt protein (231-491AA)
Immunogen Species:	Rattus norvegicus (Rat)
Buffer:	Preservative: 0.03% Proclin 300 Constituents: 50% Glycerol, 0.01M PBS, PH 7.4
Form:	Liquid

Manufacturers Statement: This final kit system is assembled and quality-released by Assay Genie Limited.

Preparation & Storage


Storage: Upon receipt, store at -20°C or -80°C. Avoid repeated freeze.
Store Bradford Reagent at Room Temperature for 1 Year.

Recommended Dilutions:	Application	Recommended Dilution
	WB	1:1000-1:5000

Protein Quantification (Optional): To quantify total protein levels, use the Bradford Reagent included in this kit. Visit <https://www.assaygenie.com/bradford-protein-assay-protocol/> to view the full protocol

Validation Data

Image

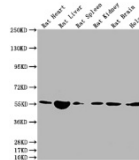

Description

Western Blot Positive WB detected in: Rat Heart tissue, Rat Liver tissue, Rat Spleen tissue, Rat Kidney tissue, Rat Brain tissue, Hela whole cell lysate All lanes: Nampt antibody at 1:2000 Secondary Goat polyclonal to rabbit IgG at 1/50000 dilution Predicted band size: 56 kDa Observed band size: 56 kDa