

Recombinant Mouse TNFSF11/RANKL/CD254 Protein

RPCB0011

Description

This high-purity Recombinant Mouse TNFSF11/RANKL/CD254 Protein is supplied as a kit for advanced applications. The kit includes Bradford Reagent to quantify total protein concentration for accurate sample normalization (Optional).

Protein Information

SKU:	RPCB0011
Contents:	10 µg, 20 µg, 50 µg, 100 µg Bradford Reagent: 1 vial (2ml)
Synonyms:	ODF, OPGL, RANKL, Ly109l, Trance, TNFSF11
Species:	Mouse
Gene ID:	21943
Expression Host:	HEK293 cells
Tags:	N-His & N-hFC
Calculated MW:	55.63 KD
Observed MW:	60-75 KD
Purification:	≥ 95 % as determined by SDS-PAGE, ≥ 95 % as determined by HPLC.
Endotoxin:	< 0.01 EU/µg of the protein by LAL method
Formulation:	Lyophilized from a 0.22 µm filtered solution of PBS, pH 7.4.
Bio- Activity:	M-CSF (50 ng/mL) and RANKL (100 ng/mL) can induced differentiation of osteoclasts from primary mouse bone marrow cells for 3 days; RANKL (100 ng/mL) can induced differentiation of osteoclasts from RAW264.7 cells for 5 days.


Manufacturers Statement - This final kit system is assembled and quality-released by Assay Genie Limited

Preparation & Storage

- Shipping:** The product is shipped at ambient temperature. Upon receipt, store it immediately at the temperature recommended below.
- Storage:** Store at -20°C. Store the lyophilized protein at -20°C to -80 °C up to 1 year from the date of receipt. After reconstitution, the protein solution is stable at -20°C for 3 months, at 2-8°C for up to 1 week.
Store Bradford Reagent at Room Temperature for 1 Year.
- Reconstitution:** Centrifuge the vial before opening. Reconstitute to a concentration of 0.1-0.5 mg/mL in sterile distilled water. Avoid vortex or vigorously pipetting the protein. For long term storage, it is recommended to add a carrier protein or stabilizer (e.g. 0.1% BSA, 5% HSA, 10% FBS or 5% Trehalose), and aliquot the reconstituted protein solution to minimize free-thaw cycles.
- Protein Quantification (Optional):** To quantify total protein levels, use the Bradford Reagent included in this kit. Visit <https://www.assaygenie.com/bradford-protein-assay-protocol/> to view the full protocol.


Validation Data

Image

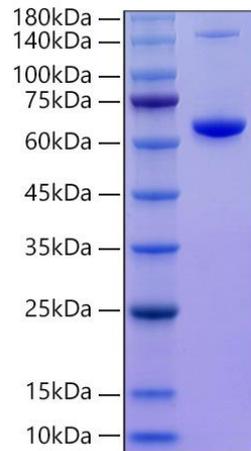

Description

Recombinant Mouse
TNFSF11/RANKL/CD254
Protein was determined by
SDS-PAGE under reducing
conditions with Coomassie
Blue.

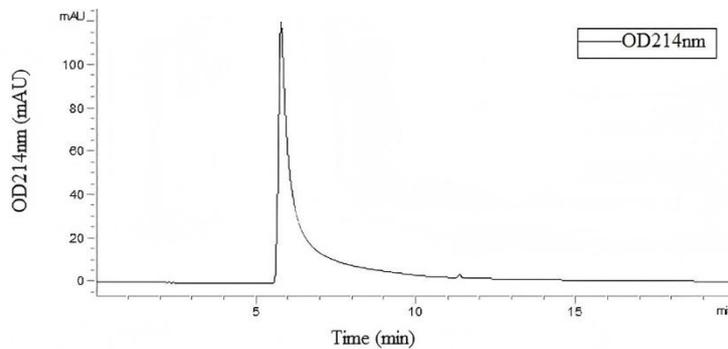

Recombinant Mouse
TNFSF11/RANKL/CD254
Protein is greater than 95% as
determined by SEC-HPLC.