

Recombinant Human IL-5 Protein

RPCB1163

Description

This high-purity Recombinant Human IL-5 Protein is supplied as a kit for advanced applications. The kit includes Bradford Reagent to quantify total protein concentration for accurate sample normalization (Optional).

Protein Information

SKU:	RPCB1163
Contents:	10 µg, 20 µg, 50 µg, 100 µg Bradford Reagent: 1 vial (2ml)
Synonyms:	IL5, EDF, IL-5, TRF
Species:	Human
Gene ID:	3567
Expression Host:	HEK293 cells
Tags:	C-His
Calculated MW:	13.99 kDa
Observed MW:	18-20 kDa
Purification:	≥ 95 % as determined by SDS-PAGE.
Endotoxin:	< 1 EU/µg of the protein by LAL method.
Formulation:	Lyophilized from a 0.22 µm filtered solution of PBS, pH 7.4.
Bio- Activity:	Measured in a cell proliferation assay using TF-1 Human erythroleukemic cells. The ED 50 for this effect is 0.5-2ng/mL.

Preparation & Storage

Shipping: The product is shipped at ambient temperature. Upon receipt, store it immediately at the temperature recommended below.


Storage: Store at -20°C. Store the lyophilized protein at -20°C to -80 °C up to 1 year from the date of receipt. After reconstitution, the protein solution is stable at -20°C for 3 months, at 2-8°C for up to 1 week.
Store Bradford Reagent at Room Temperature for 1 Year.

Reconstitution: Centrifuge the vial before opening. Reconstitute to a concentration of 0.1-0.5 mg/mL in sterile distilled water. Avoid vortex or vigorously pipetting the protein. For long term storage, it is recommended to add a carrier protein or stabilizer (e.g. 0.1% BSA, 5% HSA, 10% FBS or 5% Trehalose), and aliquot the reconstituted protein solution to minimize free-thaw cycles.

Protein Quantification (Optional): To quantify total protein levels, use the Bradford Reagent included in this kit. Visit <https://www.assaygenie.com/bradford-protein-assay-protocol/> to view the full protocol.


Validation Data

Image

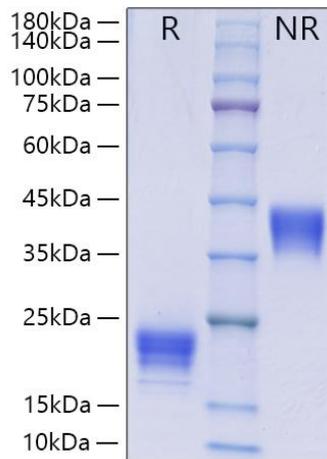

Description

Recombinant Human IL-5 Protein was determined by SDS-PAGE under reducing (R) and non-reducing (NR) conditions.

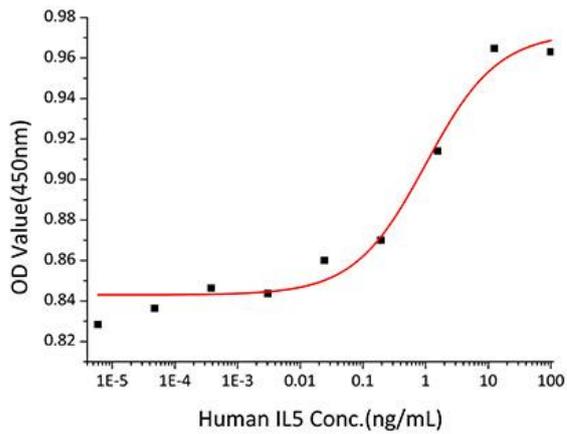

Recombinant Human IL-5 stimulates cell proliferation of the TF-1 Human erythroleukemic cells. The ED50 for this effect is 0.5-2ng/mL.